

# **Electroencephalography (EEG)**

## **Signal Processing and Analysis in Python**

**Neekesh Panchal**

November 17, 2024

## Author's Note

*For the love of education and interdisciplinary magic.*

# Contents

List of Figures	vii
I Foundations of EEG	1
1 Introduction to EEG	2
1.1 What is Electroencephalography (EEG)?	2
1.1.1 Overview of EEG Technology	2
1.1.2 Historical Development and Milestones	3
1.2 Basic Neurophysiology	3
1.2.1 Neuronal Activity and Brain Waves	3
1.2.2 Generation of EEG Signals	3
1.3 Applications of EEG	4
1.3.1 Clinical Diagnostics	4
1.3.2 Research in Neuroscience and Psychology	4
1.3.3 Brain-Computer Interfaces (BCIs)	4
1.4 Structure of the Book	4
1.4.1 Overview of Upcoming Chapters	4
1.4.2 Learning Objectives	5
2 EEG Recording Techniques and Data Structures	6
2.1 EEG Hardware Components	6
2.1.1 Electrodes and Electrode Placement Systems (10-20 System)	6
2.1.2 Amplifiers and Signal Acquisition Devices	6
2.1.3 Analog-to-Digital Conversion	8
2.2 EEG Signal Characteristics	8
2.2.1 Frequency Bands (Delta, Theta, Alpha, Beta, Gamma)	8
2.2.2 Amplitude and Phase Information	8
2.3 Data Acquisition Protocols	8
2.3.1 Setting Up an EEG Experiment	9

---

2.3.2 Subject Preparation and Safety Considerations . . . . .	9
2.4 EEG File Formats and Data Storage . . . . .	9
2.4.1 Common Formats (EDF, BDF, GDF, FIF) . . . . .	9
2.4.2 Metadata and Annotation Standards . . . . .	10
2.4.3 Data Compression and Storage Solutions . . . . .	10
II Challenges in EEG Data Analysis . . . . .	12
3 Common Problems in EEG Data . . . . .	13
3.1 Physiological Artifacts . . . . .	13
3.1.1 Eye Movements and Blinks (Electrooculogram - EOG) . . . . .	13
3.1.2 Muscle Activity (Electromyogram - EMG) . . . . .	13
3.1.3 Cardiac Signals (Electrocardiogram - ECG) . . . . .	16
3.2 Environmental and Technical Artifacts . . . . .	16
3.2.1 Electrical Noise and Interference . . . . .	16
3.2.2 Electrode Impedance Issues . . . . .	16
3.2.3 Movement Artifacts . . . . .	16
3.3 Signal Variability . . . . .	16
3.3.1 Inter-Subject and Intra-Subject Variability . . . . .	18
3.3.2 Age, Gender, and Health Factors . . . . .	18
3.4 Ethical and Privacy Concerns . . . . .	18
3.4.1 Data Anonymization . . . . .	18
3.4.2 Informed Consent and Ethical Approval . . . . .	18
III Processing EEG Data with Python . . . . .	19
4 Introduction to Python for EEG Analysis . . . . .	20
4.1 Setting Up Your Python Environment . . . . .	20
4.1.1 Python Installation (Anaconda Distribution) . . . . .	20
4.1.2 Virtual Environments and Package Management . . . . .	21
4.2 Essential Python Libraries . . . . .	21
4.2.1 NumPy and SciPy for Numerical Computing . . . . .	21
4.2.2 Matplotlib and Seaborn for Visualization . . . . .	23
4.2.3 Pandas for Data Manipulation . . . . .	23
4.2.4 MNE-Python for EEG Analysis . . . . .	23
4.3 Coding Best Practices . . . . .	23

---

4.3.1	Writing Readable and Efficient Code	23
4.3.2	Version Control with Git	24
4.3.3	Documentation and Commenting	24
<b>5</b>	<b>Reading and Understanding EEG Data in Python</b>	<b>25</b>
5.1	Introduction to MNE-Python	25
5.1.1	Overview and Capabilities	25
5.1.2	Loading Raw EEG Data	26
5.2	Importing Various EEG File Formats	26
5.2.1	EDF, BDF, GDF, and Others	26
5.2.2	Handling Data Annotations and Events	27
5.3	Exploring EEG Data Structures	27
5.3.1	Raw Objects	27
5.3.2	Epochs and Evoked Responses	27
5.3.3	Info Objects (Metadata)	28
5.4	Visualizing Raw EEG Signals	28
5.4.1	Plotting Time-Series Data	28
5.4.2	Inspecting Channel Information	29
<b>IV</b>	<b>Preprocessing and Artifact Handling</b>	<b>32</b>
<b>6</b>	<b>Data Normalization and Preprocessing</b>	<b>33</b>
6.1	Filtering Techniques	33
6.1.1	High-Pass, Low-Pass, and Band-Pass Filters	33
6.1.2	Notch Filtering for Power Line Noise	34
6.1.3	Designing FIR and IIR Filters	34
6.2	Resampling and Decimation	34
6.2.1	Changing the Sampling Rate	34
6.2.2	Anti-Aliasing Considerations	34
6.3	Baseline Correction and Detrending	35
6.3.1	Removing DC Offsets	35
6.3.2	Linear and Non-Linear Detrending Methods	35
6.4	Normalization Methods	35
6.4.1	Z-Score Normalization	35
6.4.2	Scaling Techniques	36
6.4.3	Logarithmic and Non-Linear Transformations	36

<b>7</b>	<b>Artifact Detection and Removal</b>	<b>39</b>
7.1	Identifying Artifacts . . . . .	39
7.1.1	Visual Inspection Strategies . . . . .	39
7.1.2	Automated Detection Algorithms . . . . .	40
7.2	Correction Techniques . . . . .	40
7.2.1	Independent Component Analysis (ICA) . . . . .	40
7.2.2	Regression Methods . . . . .	41
7.2.3	Adaptive Filtering . . . . .	41
7.3	Evaluating Preprocessing Effectiveness . . . . .	41
7.3.1	Quality Metrics . . . . .	43
7.3.2	Comparing Pre- and Post-Correction Data . . . . .	43
<b>V</b>	<b>Advanced EEG Signal Analysis</b>	<b>46</b>
<b>8</b>	<b>Feature Extraction and Time-Frequency Analysis</b>	<b>47</b>
8.1	Time-Domain Analysis . . . . .	47
8.1.1	Event-Related Potentials (ERPs) . . . . .	47
8.1.2	Microstate Analysis . . . . .	48
8.2	Frequency-Domain Analysis . . . . .	48
8.2.1	Power Spectral Density (PSD) . . . . .	48
8.2.2	Fast Fourier Transform (FFT) Techniques . . . . .	49
8.3	Time-Frequency Representations . . . . .	49
8.3.1	Short-Time Fourier Transform (STFT) . . . . .	49
8.3.2	Wavelet Transformations . . . . .	49
8.4	Spatial Analysis . . . . .	50
8.4.1	Source Localization . . . . .	50
8.4.2	Topographical Mapping . . . . .	50
<b>9</b>	<b>Machine Learning and EEG</b>	<b>54</b>
9.1	Introduction to Machine Learning Concepts . . . . .	54
9.1.1	Supervised vs. Unsupervised Learning . . . . .	54
9.1.2	Model Evaluation Metrics . . . . .	55
9.2	Feature Selection and Dimensionality Reduction . . . . .	55
9.2.1	Principal Component Analysis (PCA) . . . . .	55
9.2.2	t-Distributed Stochastic Neighbor Embedding (t-SNE) . . . . .	56
9.3	Classification Algorithms . . . . .	56

9.3.1 Support Vector Machines (SVM)	56
9.3.2 Random Forests	56
9.3.3 K-Nearest Neighbors (KNN)	57
9.4 Deep Learning Approaches	57
9.4.1 Convolutional Neural Networks (CNN) for EEG	57
9.4.2 Recurrent Neural Networks (RNN) and LSTM	58
9.5 Model Implementation with Python	59
9.5.1 Scikit-Learn for Traditional ML	59
9.5.2 TensorFlow and PyTorch for Deep Learning	60
<b>10 Visualizing EEG Data</b>	<b>63</b>
10.1 Time-Series Visualization	63
10.1.1 Multi-Channel Plotting	63
10.1.2 Interactive Visualization Tools	64
10.2 Spectral and Topographical Plots	64
10.2.1 Spectrograms and Heatmaps	64
10.2.2 Scalp Maps and 3D Brain Representations	65
10.3 Advanced Visualization Techniques	65
10.3.1 Connectivity Graphs	66
10.3.2 Dynamic Visualizations and Animations	66
10.4 Custom Visualization with Matplotlib and MNE-Python	67
<b>A Installation Guides and Troubleshooting</b>	<b>68</b>
A.1 Setting Up Python and MNE-Python	68
A.2 Common Issues and Solutions	68
<b>B Sample Datasets and Code Examples</b>	<b>70</b>
B.1 Accessing Public EEG Databases	70
B.2 Annotated Code Snippets for Key Procedures	70
B.2.1 Loading and Plotting Raw EEG Data	70
B.2.2 Filtering EEG Data	71
B.2.3 Epoching Data Around Events	71
B.2.4 Computing and Plotting Evoked Responses	71
<b>C Mathematical Foundations</b>	<b>72</b>
C.1 Linear Algebra and Signal Processing Basics	72
C.2 Statistical Methods in EEG Analysis	72

<b>D Additional Resources</b>	<b>73</b>
D.1 Recommended Books and Articles . . . . .	73
D.2 Online Courses and Tutorials . . . . .	73
<b>E Glossary of Terms</b>	<b>74</b>
<b>F Practice Project: Generating Synthetic EEG Data and     Implementing Clustering Algorithms</b>	<b>75</b>
F.1 Step 1: Importing Necessary Libraries . . . . .	75
F.2 Step 2: Generating Synthetic EEG Data . . . . .	76
F.3 Step 3: Applying a Bandpass Filter . . . . .	76
F.4 Step 4: Clustering the EEG Data . . . . .	77
F.5 Step 5: Visualizing the Results . . . . .	77
F.6 Step 6: Running the Program . . . . .	78



# List of Figures

2.1	International 10-20 Electrode Placement System . . . . .	7
3.1	EEG recording showing EOG artifacts due to eye blinks. . . . .	14
3.2	EEG recording displaying EMG artifacts from muscle activity. . . . .	15
3.3	EEG recording with visible ECG artifacts. . . . .	17
4.1	Anaconda Installation Process . . . . .	22
4.2	Setting Up a Virtual Environment . . . . .	22
5.1	Visualization of raw EEG signals using MNE-Python. . . . .	30
6.1	Illustration of filtering applied to EEG data, showing raw and filtered signals. . . . .	37
6.2	Comparison of Power Spectral Density of EEG signals before and after downsampling. . . . .	37
7.1	Example of artifacts in EEG data detected via visual inspection. . . . .	42
7.2	Visualization of ICA components. . . . .	42
7.3	Comparison of EEG data before and after artifact correction. . . . .	44
8.1	Power Spectral Density of EEG data, highlighting dominant frequency bands. . . . .	51
8.2	Topographical map of EEG activity at specific time points. . . . .	52
F.1	Visualized EEG signal, pre and post-filtering, and 2D PCA . . . . .	78
F.2	Running the program within the terminal . . . . .	79

# Part I

## Foundations of EEG

# Chapter 1

## Introduction to EEG

### 1.1 What is Electroencephalography (EEG)?

Electroencephalography (EEG) is a non-invasive technique for recording the brain's electrical activity. By placing electrodes on the scalp, EEG captures the collective electrical signals produced by neuronal activity, providing insights into brain function and aiding in the diagnosis of neurological disorders. This method offers high temporal resolution, making it invaluable for monitoring rapid neural dynamics.

#### 1.1.1 Overview of EEG Technology

EEG technology involves several key components:

- **Electrodes:** Sensors placed on the scalp to detect electrical activity.
- **Amplifiers:** Devices that enhance the weak electrical signals for analysis.
- **Analog-to-Digital Converters (ADCs):** Convert the analog signals into digital form for processing.
- **Recording Systems:** Software and hardware that store and display the EEG data.

The standard electrode placement system, known as the 10-20 system, ensures consistent and reproducible electrode positioning across individuals, facilitating comparative studies and clinical assessments.

### 1.1.2 Historical Development and Milestones

The origins of EEG date back to the early 20th century. In 1924, German psychiatrist Hans Berger made the first recording of human brain activity using EEG, marking a significant milestone in neuroscience. Berger's pioneering work laid the foundation for modern EEG applications, including the identification of various brain wave patterns and their association with different mental states and neurological conditions. [?]

## 1.2 Basic Neurophysiology

Understanding the neurophysiological basis of EEG is essential for interpreting its signals.

### 1.2.1 Neuronal Activity and Brain Waves

Neurons communicate through electrical impulses known as action potentials. The synchronous activity of large populations of neurons generates oscillatory patterns detectable by EEG. These patterns, termed brain waves, are categorized based on their frequency ranges:

- **Delta Waves (0.5–4 Hz):** Associated with deep sleep stages.
- **Theta Waves (4–8 Hz):** Linked to light sleep and relaxation.
- **Alpha Waves (8–13 Hz):** Present during relaxed wakefulness, especially with closed eyes.
- **Beta Waves (13–30 Hz):** Related to active thinking and concentration.
- **Gamma Waves (30–100 Hz):** Involved in higher cognitive functions.

These oscillations reflect the brain's functional states and are crucial for understanding neural processing.

### 1.2.2 Generation of EEG Signals

EEG signals primarily originate from postsynaptic potentials in cortical pyramidal neurons. When these neurons are activated synchronously, their combined electrical activity produces voltage fluctuations that propagate through the brain tissue, skull, and scalp, where they are detected by surface electrodes. The amplitude and frequency of these signals provide information about the underlying neural processes.

## 1.3 Applications of EEG

EEG has diverse applications across clinical and research domains.

### 1.3.1 Clinical Diagnostics

In clinical settings, EEG is instrumental in diagnosing and monitoring neurological disorders such as epilepsy, sleep disorders, and encephalopathies. It aids in identifying abnormal brain activity patterns, guiding treatment decisions, and evaluating therapeutic outcomes. [?]

### 1.3.2 Research in Neuroscience and Psychology

Researchers utilize EEG to study cognitive processes, sensory perception, and neural mechanisms underlying behavior. Its high temporal resolution allows for the examination of rapid neural events, making it a valuable tool in cognitive neuroscience and psychology.

### 1.3.3 Brain-Computer Interfaces (BCIs)

EEG serves as a foundation for BCIs, enabling direct communication between the brain and external devices. This technology has applications in assistive devices for individuals with motor impairments, neurorehabilitation, and even gaming.

## 1.4 Structure of the Book

This book is structured to provide a comprehensive understanding of EEG signal processing and analysis using Python.

### 1.4.1 Overview of Upcoming Chapters

- **Chapter 2:** Explores EEG recording techniques and data structures.
- **Chapter 3:** Discusses common challenges in EEG data analysis.
- **Chapter 4:** Introduces Python tools and libraries for EEG analysis.
- **Chapter 5:** Covers data normalization and preprocessing methods.
- **Chapter 6:** Focuses on artifact detection and removal techniques.
- **Chapter 7:** Delves into advanced EEG signal analysis methods.

- **Chapter 8:** Examines machine learning applications in EEG analysis.
- **Chapter 9:** Provides guidance on visualizing EEG data.
- **Chapter 10:** Presents real-world applications and case studies.
- **Chapter 11:** Discusses emerging trends and ethical considerations in EEG research.

### 1.4.2 Learning Objectives

By the end of this book, readers will be able to:

- Understand the principles of EEG recording and data acquisition.
- Identify and address common issues in EEG data.
- Utilize Python for effective EEG data processing and analysis.
- Apply advanced signal processing techniques to extract meaningful information.
- Implement machine learning models for EEG data classification and prediction.
- Visualize EEG data to interpret and present findings effectively.
- Recognize ethical considerations and emerging trends in EEG research.

This structured approach aims to equip readers with both theoretical knowledge and practical skills in EEG signal processing and analysis.

# Chapter 2

## EEG Recording Techniques and Data Structures

### 2.1 EEG Hardware Components

Electroencephalography (EEG) relies on specialized hardware to capture and record the brain's electrical activity. This section explores the primary components involved in EEG data acquisition.

#### 2.1.1 Electrodes and Electrode Placement Systems (10-20 System)

EEG electrodes are sensors placed on the scalp to detect electrical signals generated by neuronal activity. The International 10-20 system is the standard method for electrode placement, ensuring consistency and reproducibility across recordings. In this system, electrodes are positioned based on the relative distances between anatomical landmarks on the skull, with placements at 10% or 20% intervals. This standardized approach facilitates accurate localization of brain activity and comparison across studies.

#### 2.1.2 Amplifiers and Signal Acquisition Devices

The electrical signals detected by EEG electrodes are typically in the microvolt range and susceptible to noise. Amplifiers are essential for enhancing these weak signals to levels suitable for analysis. High-quality EEG amplifiers possess characteristics such as high input impedance, low noise, and a wide dynamic range to preserve signal integrity. Modern signal acquisition devices often integrate amplifiers with analog-to-digital converters (ADCs) and interfaces for data storage and visualization.

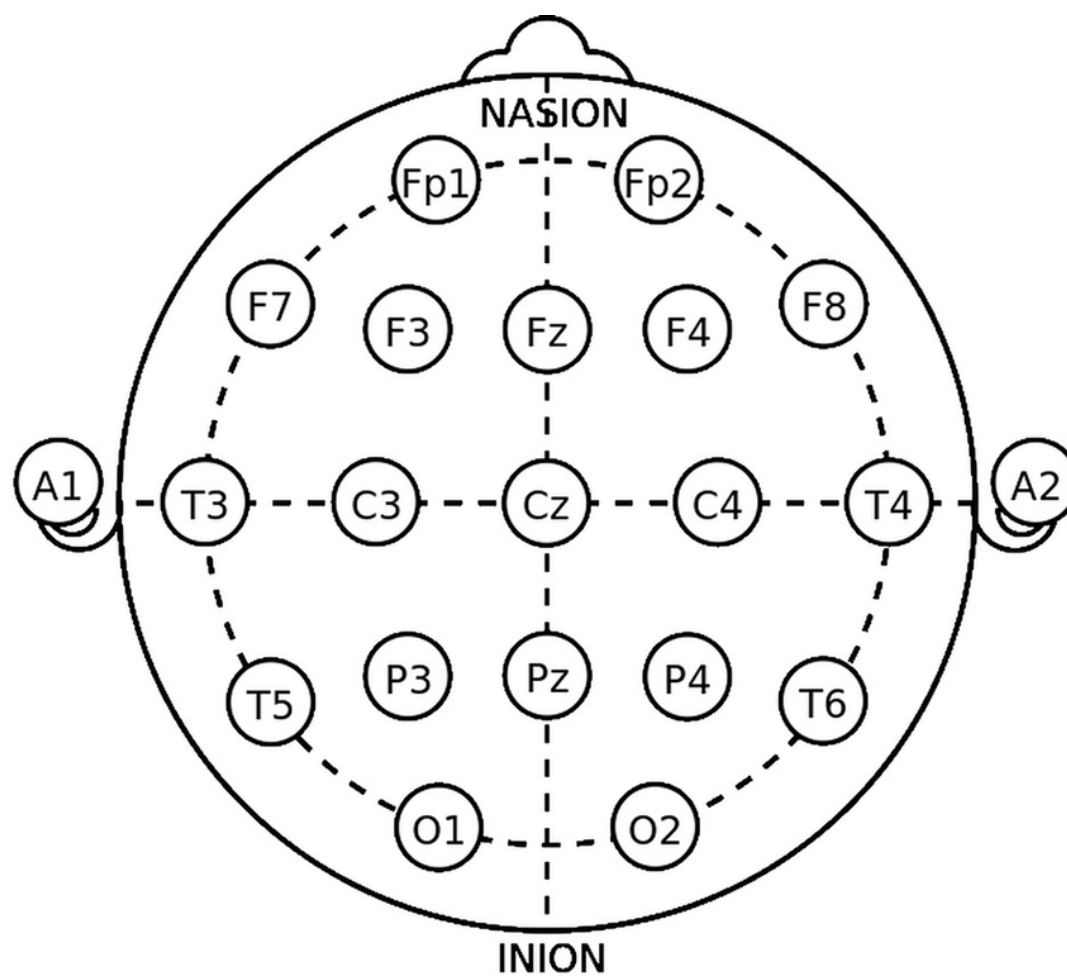


Figure 2.1: International 10-20 Electrode Placement System



### 2.1.3 Analog-to-Digital Conversion

Analog-to-digital converters (ADCs) transform the continuous analog EEG signals into discrete digital data for processing and analysis. Key parameters of ADCs include sampling rate and resolution. The sampling rate must be sufficiently high to capture the fastest components of the EEG signal, adhering to the Nyquist theorem. A higher resolution allows for more precise representation of signal amplitude, enhancing the quality of the recorded data.

## 2.2 EEG Signal Characteristics

Understanding the fundamental characteristics of EEG signals is crucial for effective analysis and interpretation.

### 2.2.1 Frequency Bands (Delta, Theta, Alpha, Beta, Gamma)

EEG signals are composed of oscillations across various frequency ranges, each associated with different cognitive and physiological states:

- **Delta (0.5–4 Hz):** Predominant during deep sleep stages.
- **Theta (4–8 Hz):** Linked to light sleep, relaxation, and meditative states.
- **Alpha (8–13 Hz):** Observed during relaxed wakefulness, especially with closed eyes.
- **Beta (13–30 Hz):** Associated with active thinking, concentration, and alertness.
- **Gamma (30–100 Hz):** Related to higher cognitive functions, including perception and consciousness.

### 2.2.2 Amplitude and Phase Information

The amplitude of EEG signals reflects the power or intensity of brain activity within specific frequency bands, while phase information pertains to the timing relationships between oscillatory components. Analyzing both amplitude and phase provides insights into neural synchrony and connectivity, offering a comprehensive understanding of brain dynamics.

## 2.3 Data Acquisition Protocols

Implementing standardized protocols during EEG data acquisition ensures the reliability and validity of the recordings.

### 2.3.1 Setting Up an EEG Experiment

Establishing a controlled environment is essential for minimizing artifacts and obtaining high-quality EEG data. This involves:

- **Environmental Control:** Utilizing a quiet, electrically shielded room to reduce external noise.
- **Equipment Calibration:** Regularly calibrating EEG equipment to maintain accuracy.
- **Protocol Standardization:** Following consistent procedures for electrode application, impedance checking, and data recording.

### 2.3.2 Subject Preparation and Safety Considerations

Proper preparation of subjects is vital for both data quality and participant safety:

- **Skin Preparation:** Cleaning the scalp to reduce impedance and enhance signal quality.
- **Electrode Application:** Ensuring secure and comfortable electrode placement to prevent movement artifacts.
- **Safety Measures:** Monitoring for any adverse reactions and ensuring that all equipment is properly grounded to prevent electrical hazards.

## 2.4 EEG File Formats and Data Storage

Efficient storage and management of EEG data are facilitated by standardized file formats and data structures.

### 2.4.1 Common Formats (EDF, BDF, GDF, FIF)

Several file formats are commonly used for storing EEG data:

- **European Data Format (EDF):** A widely adopted standard for storing bioelectrical signals, known for its simplicity and compatibility.
- **BioSemi Data Format (BDF):** Similar to EDF but supports higher resolution data, commonly used with BioSemi systems.
- **General Data Format (GDF):** An extension of EDF, incorporating additional metadata and event information.

- **Functional Imaging Format (FIF):** Used by the MNE software suite, capable of storing EEG and MEG data along with associated metadata.

### 2.4.2 Metadata and Annotation Standards

Incorporating metadata and annotations enhances the utility of EEG datasets:

- **Metadata:** Includes information about the recording conditions, subject demographics, and equipment settings.
- **Annotations:** Markers indicating events, artifacts, or specific time points of interest within the data.

Adhering to standards such as the Brain Imaging Data Structure (BIDS) for EEG promotes consistency and facilitates data sharing and collaboration.

### 2.4.3 Data Compression and Storage Solutions

EEG datasets can be large, necessitating efficient storage solutions:

- **Compression Techniques:** Methods such as lossless compression reduce file size without sacrificing data integrity.
- **Storage Solutions:** Utilizing robust data management systems and cloud storage options ensures data security and accessibility.

Implementing these strategies aids in managing storage requirements and facilitates efficient data retrieval for analysis.

# Bibliography

- [1] Jurcak, V., Tsuzuki, D., & Dan, I. (2007). 10/20, 10/10, and 10/5 systems revisited: Their validity as relative head-surface-based positioning systems. *NeuroImage*, 34(4), 1600-1611.
- [2] Kemp, B., Värri, A., Rosa, A. C., Nielsen, K. D., & Gade, J. (1992). A simple format for exchange of digitized polygraphic recordings. *Electroencephalography and Clinical Neurophysiology*, 82(5), 391-393.
- [3] Gorgolewski, K. J., Auer, T., Calhoun, V. [::contentReference\[oaicite:0\]index=0](#)

## **Part II**

# **Challenges in EEG Data Analysis**

# Chapter 3

## Common Problems in EEG Data

Electroencephalography (EEG) is a powerful tool for monitoring brain activity, but its effectiveness can be compromised by various artifacts and challenges. This chapter delves into the common issues encountered in EEG data analysis, including physiological and environmental artifacts, signal variability, and ethical considerations.

### 3.1 Physiological Artifacts

Physiological artifacts originate from the subject's own biological activities, which can interfere with the EEG signals.

#### 3.1.1 Eye Movements and Blinks (Electrooculogram - EOG)

Eye movements and blinks generate electrical potentials that can overshadow the brain's electrical activity. These artifacts are typically observed as large, slow waves in the EEG, predominantly affecting the frontal electrodes. The corneo-retinal dipole model explains this phenomenon, where the cornea is positively charged relative to the retina, and eye movements alter the electric field detected by scalp electrodes.

#### 3.1.2 Muscle Activity (Electromyogram - EMG)

Muscle contractions, especially from facial and neck muscles, produce high-frequency artifacts that can contaminate EEG recordings. These artifacts are characterized by increased power in the gamma frequency band (30–100 Hz) and can obscure neural signals.

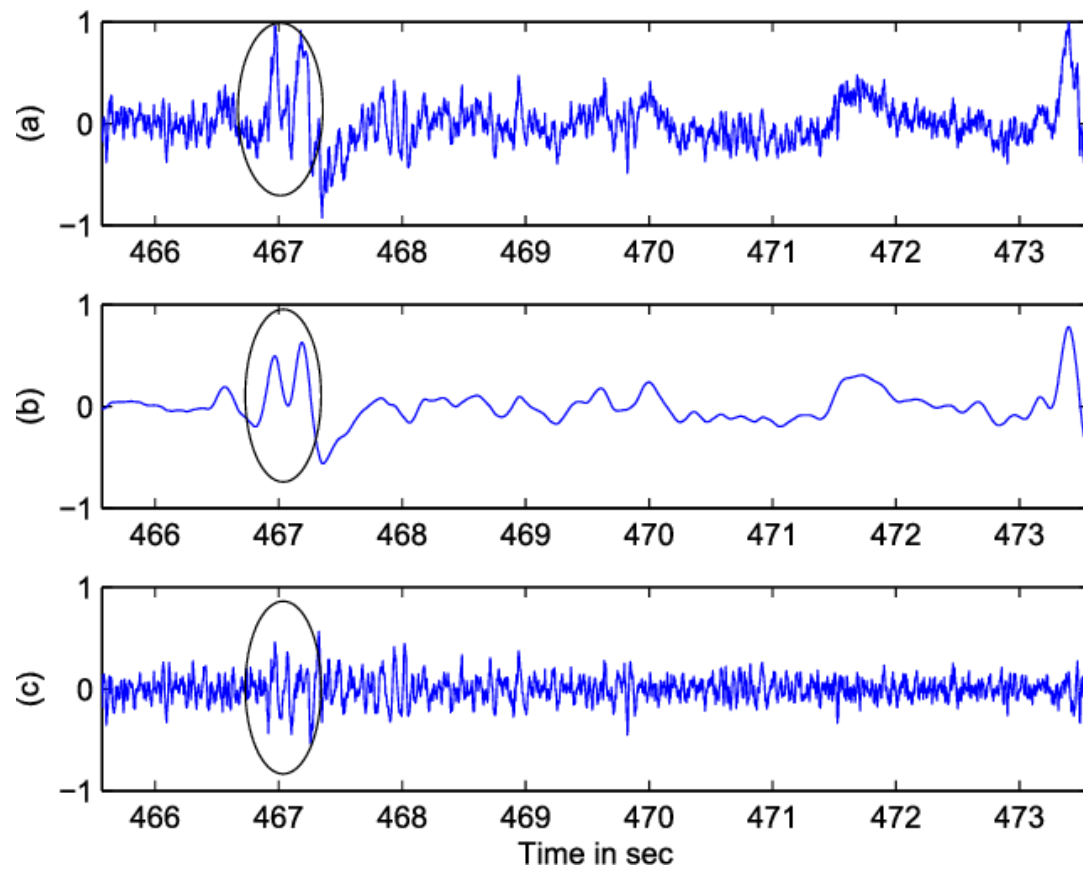


Figure 3.1: EEG recording showing EOG artifacts due to eye blinks.

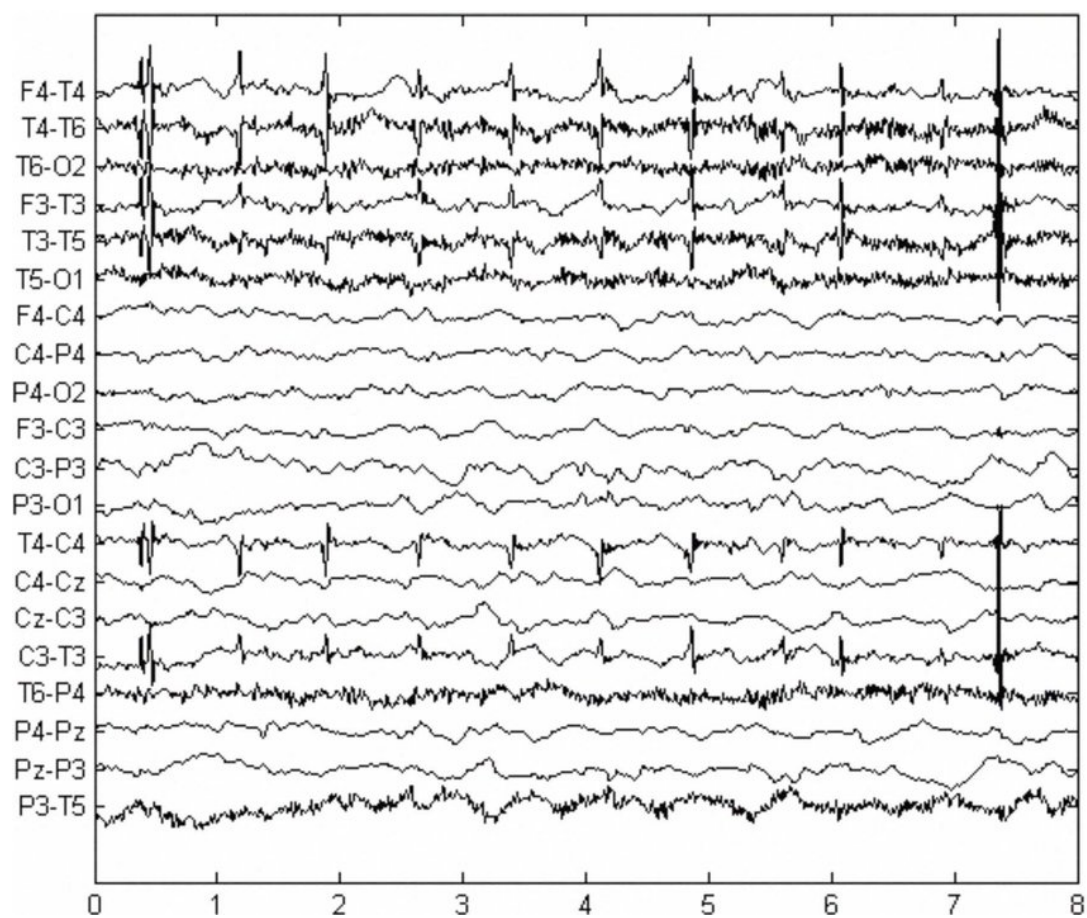


Figure 3.2: EEG recording displaying EMG artifacts from muscle activity.



### **3.1.3 Cardiac Signals (Electrocardiogram - ECG)**

The heart's electrical activity can introduce artifacts into EEG recordings, known as electrocardiographic (ECG) artifacts. These are typically rhythmic and can be mistaken for brain activity if not properly identified and removed.

## **3.2 Environmental and Technical Artifacts**

External factors and equipment-related issues can also introduce artifacts into EEG data.

### **3.2.1 Electrical Noise and Interference**

Electrical devices and power lines can emit electromagnetic fields that interfere with EEG recordings, resulting in artifacts, often observed as 50 or 60 Hz oscillations, depending on the local power supply frequency. Proper grounding and shielding of EEG equipment are essential to minimize these artifacts.

### **3.2.2 Electrode Impedance Issues**

High or uneven electrode impedance can lead to poor signal quality and increased susceptibility to noise. Regular impedance checks and proper electrode application techniques are crucial to ensure reliable EEG recordings.

### **3.2.3 Movement Artifacts**

Subject movements, such as head shifts or body movements, can cause artifacts in EEG data. These artifacts are typically low-frequency signals that can mask the underlying brain activity.

## **3.3 Signal Variability**

EEG signals can vary due to individual differences and external factors.



Figure 3.3: EEG recording with visible ECG artifacts.

### **3.3.1 Inter-Subject and Intra-Subject Variability**

Variations between different individuals (inter-subject) and within the same individual over time (intra-subject) can affect EEG data. Factors such as anatomical differences, mental state, and fatigue levels contribute to this variability, posing challenges for data analysis and interpretation.

### **3.3.2 Age, Gender, and Health Factors**

Demographic factors like age and gender, as well as health conditions, can influence EEG signals. For instance, aging is associated with changes in EEG rhythms, and neurological disorders can alter typical EEG patterns.

## **3.4 Ethical and Privacy Concerns**

Handling EEG data involves ethical considerations to protect participants' rights and privacy.

### **3.4.1 Data Anonymization**

Ensuring that EEG data cannot be traced back to individual participants is crucial for privacy. Techniques such as removing personal identifiers and using anonymized codes are standard practices.

### **3.4.2 Informed Consent and Ethical Approval**

Obtaining informed consent from participants and securing ethical approval from relevant committees are fundamental steps in conducting EEG research. These measures ensure that participants are aware of the study's purpose and their rights, and that the research adheres to ethical standards.

Addressing these challenges is essential for obtaining high-quality EEG data and ensuring the validity and reliability of subsequent analyses.

## Part III

# Processing EEG Data with Python

# Chapter 4

## Introduction to Python for EEG Analysis

### 4.1 Setting Up Your Python Environment

To effectively analyze EEG data using Python, it's essential to establish a robust and organized development environment. This section guides you through the installation of Python, the creation of virtual environments, and the management of necessary packages.

#### 4.1.1 Python Installation (Anaconda Distribution)

The Anaconda distribution is a comprehensive platform that simplifies Python installation and package management, making it particularly suitable for scientific computing and data analysis tasks.

- **Download Anaconda:** Visit the official Anaconda website at <https://www.anaconda.com/products/distribution> and download the installer compatible with your operating system.
- **Install Anaconda:** Run the downloaded installer and follow the on-screen instructions to complete the installation process.
- **Verify Installation:** Open a terminal or command prompt and type `conda -version` to confirm that Anaconda is installed correctly.

### 4.1.2 Virtual Environments and Package Management

Managing dependencies is crucial in Python projects to ensure compatibility and reproducibility. Virtual environments allow you to isolate project-specific packages, preventing conflicts between different projects.

- **Create a Virtual Environment:** Use the following command to create a new virtual environment named `eeg_env`:

```
conda create --name eeg_env python=3.8
```

- **Activate the Environment:** Activate the newly created environment with:

```
conda activate eeg_env
```

- **Install Necessary Packages:** Within the activated environment, install required packages using `pip` or `conda`. For example:

```
pip install numpy scipy matplotlib mne
```

- **Deactivate the Environment:** To exit the virtual environment, use:

```
conda deactivate
```

## 4.2 Essential Python Libraries

Several Python libraries are indispensable for EEG data analysis. This section introduces key libraries and their primary functionalities.

### 4.2.1 NumPy and SciPy for Numerical Computing

- **NumPy:** Provides support for large, multi-dimensional arrays and matrices, along with a collection of mathematical functions to operate on these arrays. It's fundamental for numerical computations in Python.

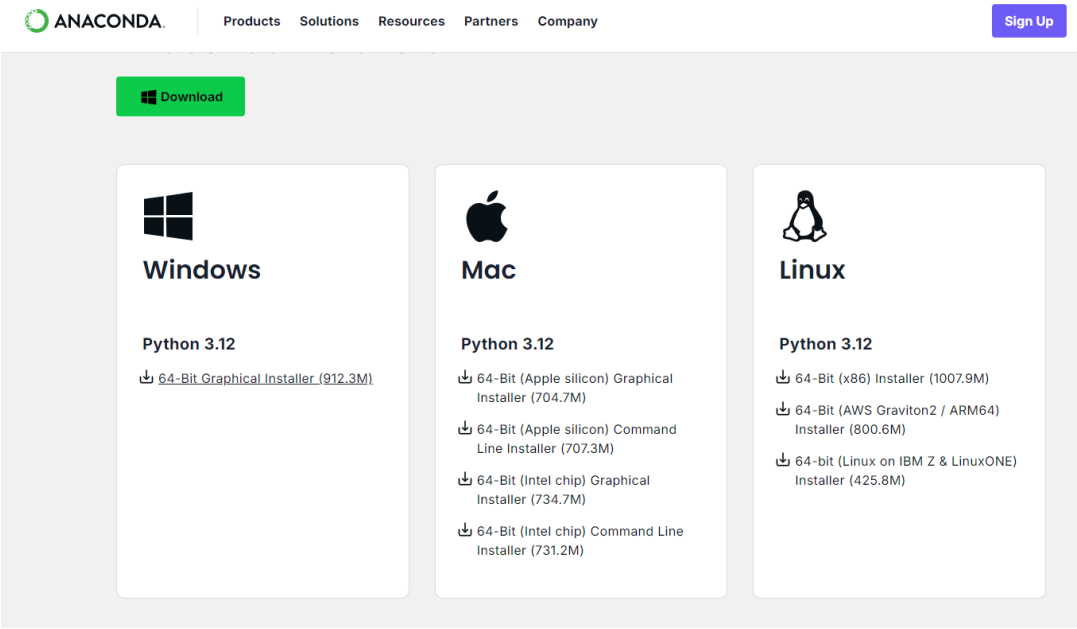


Figure 4.1: Anaconda Installation Process

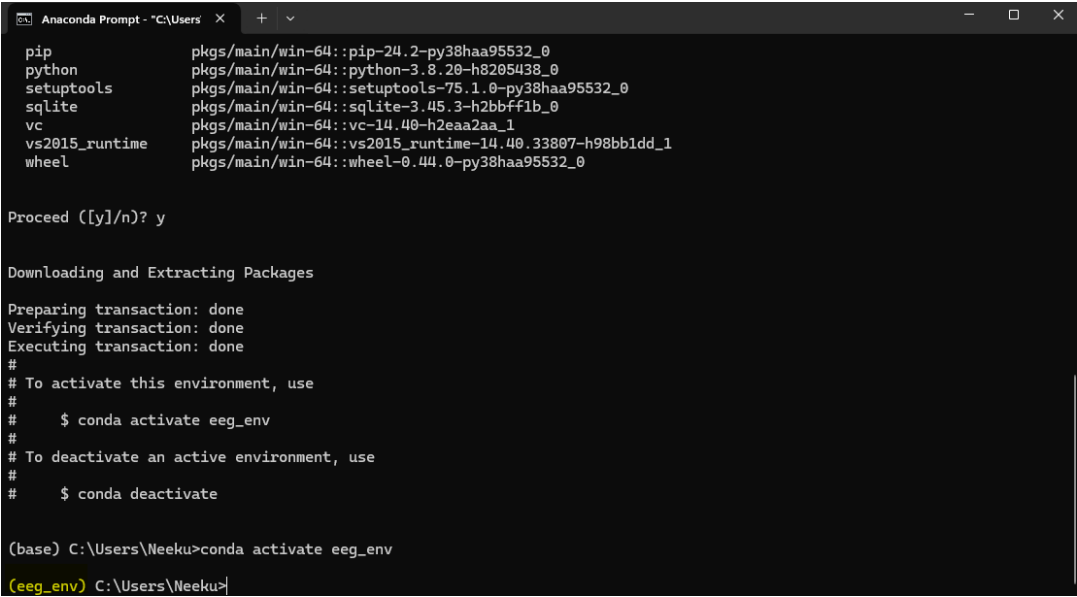


Figure 4.2: Setting Up a Virtual Environment

- **SciPy:** Builds on NumPy by adding a collection of algorithms and high-level commands for data manipulation and analysis, including signal processing functions essential for EEG data analysis.

### 4.2.2 Matplotlib and Seaborn for Visualization

- **Matplotlib:** A versatile plotting library that enables the creation of static, interactive, and animated visualizations in Python. It's widely used for generating plots of EEG signals.
- **Seaborn:** Built on top of Matplotlib, Seaborn provides a high-level interface for drawing attractive and informative statistical graphics, facilitating the visualization of complex data relationships.

### 4.2.3 Pandas for Data Manipulation

- **Pandas:** Offers data structures and functions needed to work with structured data seamlessly. It provides tools for reading and writing data between in-memory data structures and different formats, including CSV and Excel.

### 4.2.4 MNE-Python for EEG Analysis

- **MNE-Python:** A comprehensive library specifically designed for processing and analyzing EEG and MEG data. It includes tools for data preprocessing, visualization, source localization, and statistical analysis, making it a cornerstone for EEG research in Python.

## 4.3 Coding Best Practices

Adhering to coding best practices enhances code readability, maintainability, and collaboration. This section outlines key practices to follow in your EEG analysis projects.

### 4.3.1 Writing Readable and Efficient Code

- **Use Descriptive Variable Names:** Choose meaningful names that convey the purpose of the variable.



- **Follow PEP 8 Guidelines:** Adhere to Python's style guide to maintain consistency in your codebase.
- **Modularize Your Code:** Break down your code into functions and classes to promote reusability and clarity.

### 4.3.2 Version Control with Git

- **Initialize a Git Repository:** Use `git init` to start version control in your project directory.
- **Commit Changes Regularly:** Save snapshots of your code at logical points to track progress and facilitate rollback if necessary.
- **Use Branches for Features:** Develop new features or experiments in separate branches to keep the main codebase stable.

### 4.3.3 Documentation and Commenting

- **Docstrings:** Include docstrings in your functions and classes to describe their purpose and usage.
- **Inline Comments:** Add comments to explain complex or non-obvious parts of your code.
- **Maintain a README File:** Provide an overview of your project, including setup instructions and usage examples.

By setting up a well-structured Python environment and adhering to best coding practices, you lay a solid foundation for effective and efficient EEG data analysis.

# Chapter 5

## Reading and Understanding EEG Data in Python

### 5.1 Introduction to MNE-Python

Electroencephalography (EEG) data analysis necessitates robust tools capable of handling complex datasets. MNE-Python is an open-source Python package specifically designed for processing, analyzing, and visualizing EEG and magnetoencephalography (MEG) data. It offers a comprehensive suite of functions that facilitate tasks ranging from data importation to advanced signal processing and source localization [1].

#### 5.1.1 Overview and Capabilities

MNE-Python provides a wide array of features tailored for EEG data analysis:

- **Data Importation:** Supports various EEG data formats, including EDF, BDF, GDF, and FIF, enabling seamless integration with different recording systems [1].
- **Preprocessing Tools:** Offers functions for filtering, resampling, artifact detection, and correction to enhance data quality [4].
- **Data Structures:** Utilizes specialized objects such as Raw, Epochs, and Evoked to represent continuous, segmented, and averaged data, respectively [1].
- **Visualization Utilities:** Includes capabilities for plotting time-series data, topographical maps, and 3D brain representations to facilitate data interpretation [1].

- **Advanced Analysis Techniques:** Supports time-frequency analysis, source localization, and connectivity analysis, catering to complex research needs [3].

These features make MNE-Python a versatile and powerful tool for EEG data analysis.

### 5.1.2 Loading Raw EEG Data

To begin analyzing EEG data with MNE-Python, the initial step involves loading the raw data into the environment. This is accomplished using the `mne.io.read_raw_*` functions, which are tailored to handle specific file formats. For instance, to load a raw EDF file:

```
import mne
raw = mne.io.read_raw_edf('path_to_file.edf', preload=True)
```

The `preload=True` parameter loads the data into memory, allowing for efficient processing. Once loaded, the `raw` object contains the continuous EEG data along with associated metadata, serving as the foundation for subsequent analysis steps.

## 5.2 Importing Various EEG File Formats

EEG data is stored in multiple file formats, each with unique characteristics. MNE-Python supports a broad spectrum of these formats, facilitating flexibility in data handling.

### 5.2.1 EDF, BDF, GDF, and Others

The following table summarizes the functions used to read different EEG file formats in MNE-Python:

File Format	MNE-Python Function
EDF (European Data Format)	<code>mne.io.read_raw_edf</code>
BDF (BioSemi Data Format)	<code>mne.io.read_raw_bdf</code>
GDF (General Data Format)	<code>mne.io.read_raw_gdf</code>
FIF (Functional Imaging Format)	<code>mne.io.read_raw_fif</code>

For example, to load a BDF file:

```
raw = mne.io.read_raw_bdf('path_to_file.bdf', preload=True)
```

This versatility ensures compatibility with various EEG recording systems and facilitates seamless data integration [3].

### 5.2.2 Handling Data Annotations and Events

Annotations in EEG data mark specific events or periods of interest, such as stimulus presentations or artifacts. MNE-Python provides tools to manage these annotations effectively.

To access existing annotations in a raw data object:

```
annotations = raw.annotations
```

To add a new annotation:

```
raw.annotations.append(onset=10.0, duration=2.0, description='Blink')
```

This approach allows for precise marking and handling of events within the EEG data, which is crucial for accurate analysis [4].

## 5.3 Exploring EEG Data Structures

MNE-Python employs specialized data structures to represent different forms of EEG data, each serving a specific purpose in the analysis pipeline.

### 5.3.1 Raw Objects

The `Raw` object encapsulates continuous EEG data, including all channels and time points. It provides methods for data inspection, filtering, and plotting. For instance, to plot the raw data:

```
raw.plot()
```

This visualization aids in identifying artifacts and assessing data quality before proceeding with further analysis [3].

### 5.3.2 Epochs and Evoked Responses

EEG data is often segmented into epochs—time windows around specific events—to analyze event-related potentials (ERPs). In MNE-Python, the `Epochs` object represents these segments.

To create epochs:

```
events = mne.find_events(raw)
epochs = mne.Epochs(raw, events, event_id=1, tmin=-0.2, tmax=0.5, baseline=(None,
```

Averaging these epochs yields an **Evoked** object, representing the ERP:

```
evoked = epochs.average()
evoked.plot()
```

This process facilitates the examination of neural responses to specific stimuli or events [1].

### 5.3.3 Info Objects (Metadata)

The **Info** object in MNE-Python stores metadata about the EEG recording, including channel names, types, sampling frequency, and measurement information. Accessing this information is straightforward:

```
info = raw.info
print(info)
```

Understanding the metadata is essential for accurate data interpretation and ensures that subsequent analyses are appropriately configured [3].

## 5.4 Visualizing Raw EEG Signals

Effective visualization is crucial for interpreting EEG data and identifying potential issues. MNE-Python offers robust plotting functions to facilitate this process.

### 5.4.1 Plotting Time-Series Data

To visualize the time-series data of the raw EEG signals:

```
raw.plot()
```

This function provides an interactive plot where channels can be selected or deselected, and specific time windows can be examined in detail [1].

### 5.4.2 Inspecting Channel Information

To inspect the information of individual channels:

```
print(raw.info[ 'ch_names' ])
```

This command lists all channel names, allowing for verification and selection of channels for further analysis.

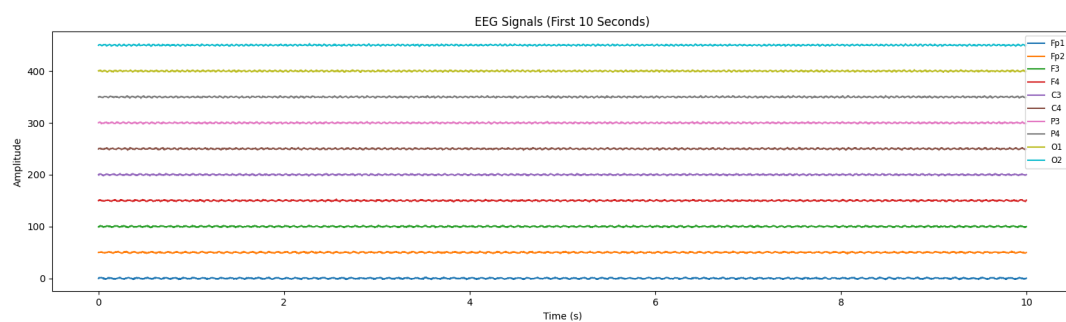


Figure 5.1: Visualization of raw EEG signals using MNE-Python.

# Bibliography

- [1] Gramfort, A., et al. (2013). MEG and EEG data analysis with MNE-Python. *Frontiers in Neuroscience*, 7, 267.
- [2] MNE-Python Documentation. (2023). Retrieved from <https://mne.tools/stable/index.html>.
- [3] Cohen, M. X. (2014). *Analyzing Neural Time Series Data: Theory and Practice*. MIT Press.
- [4] Delorme, A., Makeig, S. (2004). EEGLAB: An open-source toolbox for analysis of single-trial EEG dynamics. *Journal of Neuroscience Methods*, 134(1), 9–21.



## Part IV

# Preprocessing and Artifact Handling

# Chapter 6

## Data Normalization and Preprocessing

Data preprocessing is a critical step in EEG analysis, ensuring that the data is clean, interpretable, and ready for further processing. This chapter delves into various preprocessing techniques, including filtering, resampling, baseline correction, and normalization methods.

### 6.1 Filtering Techniques

Filtering is essential for removing unwanted noise and isolating specific frequency bands of interest in EEG data. It involves applying mathematical operations to the signal to suppress or enhance certain frequencies.

#### 6.1.1 High-Pass, Low-Pass, and Band-Pass Filters

- **High-Pass Filter:** Removes low-frequency noise, such as drift or baseline wander, by allowing frequencies above a specified cutoff to pass. Commonly used for frequencies above 0.5 Hz.
- **Low-Pass Filter:** Eliminates high-frequency noise, such as muscle artifacts, by allowing frequencies below a specified cutoff to pass. Often used with a cutoff around 30–50 Hz.
- **Band-Pass Filter:** Combines high-pass and low-pass filters to isolate a specific range of frequencies, such as the alpha band (8–13 Hz).

Filtering can be implemented in MNE-Python as follows:

```
# Apply a band-pass filter between 1 and 30 Hz  
raw.filter(l_freq=1.0, h_freq=30.0)
```

### 6.1.2 Notch Filtering for Power Line Noise

Power line noise (50/60 Hz) is a common artifact in EEG data. A notch filter effectively removes this interference by suppressing frequencies around the power line frequency.

```
# Apply a notch filter at 50 Hz  
raw.notch_filter(freqs=50.0)
```

### 6.1.3 Designing FIR and IIR Filters

- **Finite Impulse Response (FIR) Filters:** Non-recursive filters that are stable and do not introduce phase distortions, making them suitable for EEG preprocessing.
- **Infinite Impulse Response (IIR) Filters:** Recursive filters that are computationally efficient but may introduce phase distortions.

MNE-Python uses FIR filters by default, but users can switch to IIR filters if necessary:

```
# Configure and apply an IIR filter  
iir_params = mne.filter.create_filter(raw.get_data(), raw.info['sfreq'],  
                                     l_freq=1.0, h_freq=30.0, method='iir')
```

## 6.2 Resampling and Decimation

Resampling involves changing the sampling rate of the EEG data, often to reduce the data size or to match the requirements of specific analysis methods.

### 6.2.1 Changing the Sampling Rate

The sampling rate can be adjusted using the `resample` method:

```
# Resample data to 100 Hz  
raw.resample(sfreq=100)
```

### 6.2.2 Anti-Aliasing Considerations

Before resampling, an anti-aliasing filter is applied to prevent high-frequency components from aliasing into lower frequencies. MNE-Python handles this automatically during resampling.

## 6.3 Baseline Correction and Detrending

Baseline correction and detrending are essential for removing systematic biases and trends in EEG data.

### 6.3.1 Removing DC Offsets

A DC offset is a constant shift in the signal baseline, often caused by equipment artifacts. Removing this offset ensures that the data centers around zero.

```
# Apply baseline correction  
epochs.apply_baseline(baseline=(None, 0))
```

### 6.3.2 Linear and Non-Linear Detrending Methods

Detrending removes linear or non-linear trends from the data:

- **Linear Detrending:** Removes linear trends using a least-squares fit.
- **Non-Linear Detrending:** Removes trends using polynomial or spline fitting.

```
# Apply linear detrending  
raw._data = mne.preprocessing.compute_proj(raw)
```

## 6.4 Normalization Methods

Normalization standardizes EEG data, facilitating comparison across channels or subjects.

### 6.4.1 Z-Score Normalization

Z-score normalization scales data based on its mean and standard deviation, ensuring that it has a mean of 0 and a standard deviation of 1:

```
# Apply Z-score normalization  
data = raw.get_data()  
normalized_data = (data - np.mean(data, axis=1, keepdims=True)) / np.std(data, ax
```

### 6.4.2 Scaling Techniques

Scaling involves adjusting data to a specific range, such as  $[0, 1]$ :

```
# Apply Min-Max scaling  
from sklearn.preprocessing import MinMaxScaler  
scaler = MinMaxScaler()  
scaled_data = scaler.fit_transform(data.T).T
```

### 6.4.3 Logarithmic and Non-Linear Transformations

Logarithmic or non-linear transformations reduce the impact of large amplitudes:

```
# Apply logarithmic transformation  
log_data = np.log1p(data)
```

These preprocessing steps are critical for ensuring that EEG data is clean, interpretable, and ready for subsequent analysis.

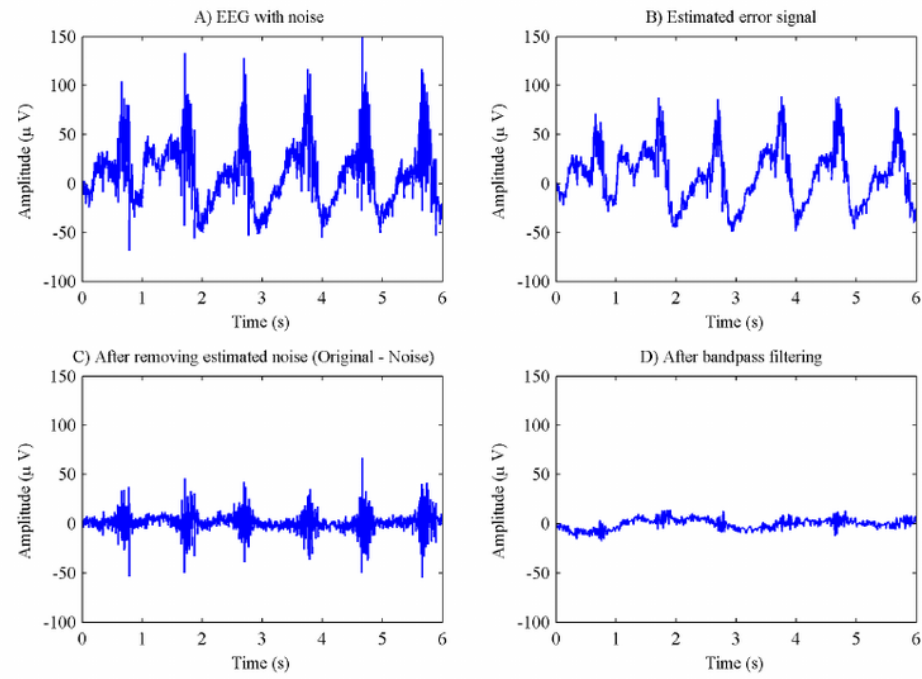


Figure 6.1: Illustration of filtering applied to EEG data, showing raw and filtered signals.

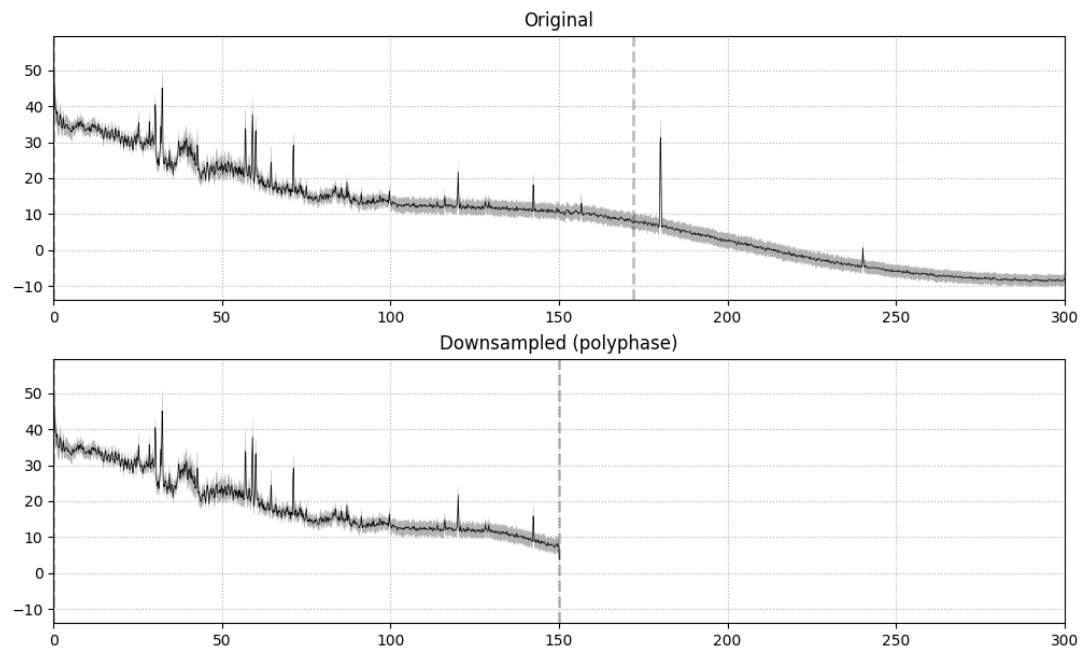


Figure 6.2: Comparison of Power Spectral Density of EEG signals before and after downsampling.

# Bibliography

- [1] Gramfort, A., et al. (2013). MEG and EEG data analysis with MNE-Python. *Frontiers in Neuroscience*, 7, 267.
- [2] MNE-Python Documentation. (2023). Retrieved from <https://mne.tools/stable/index.html>.
- [3] Cohen, M. X. (2014). *Analyzing Neural Time Series Data: Theory and Practice*. MIT Press.
- [4] Delorme, A., Makeig, S. (2004). EEGLAB: An open-source toolbox for analysis of single-trial EEG dynamics. *Journal of Neuroscience Methods*, 134(1), 9–21.

# Chapter 7

## Artifact Detection and Removal

EEG data is often contaminated by various artifacts, which can obscure meaningful neural signals. This chapter explores methods for identifying and correcting artifacts, including both visual and automated strategies, and evaluates preprocessing effectiveness.

### 7.1 Identifying Artifacts

Artifacts in EEG data can originate from physiological sources, such as eye movements or muscle activity, or from external sources, such as electrical interference.

#### 7.1.1 Visual Inspection Strategies

Visual inspection involves manually reviewing EEG data to identify artifacts based on their characteristic patterns:

- **Eye Blinks and Movements:** Large, slow waves typically affecting frontal channels.
- **Muscle Artifacts:** High-frequency noise due to muscle contractions, especially in temporal regions.
- **Line Noise:** Regular oscillations at 50 or 60 Hz caused by electrical interference.

This process is often performed using interactive plotting tools provided by MNE-Python:

```
# Plot raw data for visual inspection  
raw.plot()
```



### 7.1.2 Automated Detection Algorithms

Automated algorithms provide an efficient and objective means of identifying artifacts. Common methods include:

- **Amplitude Thresholding:** Flags segments where signal amplitude exceeds a predefined threshold.
- **Statistical Measures:** Detects artifacts based on deviations from the mean or standard deviation.
- **Machine Learning Approaches:** Classifies artifacts using supervised or unsupervised learning.

In MNE-Python, annotations can be generated automatically:

```
# Automatically detect bad segments
annotations = mne.preprocessing.annotate_raw(raw, flat=1e-6, duration=0.5)
raw.set_annotations(annotations)
```

## 7.2 Correction Techniques

Artifact correction involves modifying the EEG data to reduce or eliminate the influence of artifacts.

### 7.2.1 Independent Component Analysis (ICA)

**Theory and Implementation** ICA decomposes multichannel EEG data into statistically independent components. Artifacts, such as eye blinks, are often isolated into specific components, which can then be removed without affecting the underlying brain activity.

**Applying ICA with MNE-Python** MNE-Python provides robust ICA tools for artifact removal:

```
from mne.preprocessing import ICA

# Fit ICA on the raw data
ica = ICA(n_components=20, random_state=97)
ica.fit(raw)
```

```
# Plot ICA components for manual selection
ica.plot_components()

# Remove components associated with artifacts
ica.exclude = [0, 1] # Example: Exclude components 0 and 1
raw_ica = ica.apply(raw)
```

### 7.2.2 Regression Methods

Regression methods use reference channels (e.g., EOG for eye artifacts) to model and subtract the artifact from the EEG data. This approach assumes a linear relationship between the artifact and the recorded signal.

```
# Create a regression model for EOG artifact correction
from mne.preprocessing import create_eog_epochs

eog_epochs = create_eog_epochs(raw) # Extract EOG epochs
eog_average = eog_epochs.average() # Average the EOG artifact
raw_corrected = raw.subtract_evoked(eog_average)
```

### 7.2.3 Adaptive Filtering

Adaptive filters dynamically adjust their parameters to minimize the effect of artifacts. These filters are particularly useful for line noise and other repetitive artifacts.

```
# Apply adaptive filtering (e.g., notch filter for line noise)
raw.notch_filter(freqs=50.0)
```

## 7.3 Evaluating Preprocessing Effectiveness

Assessing the effectiveness of preprocessing is critical for ensuring that artifact removal does not distort the neural signals.

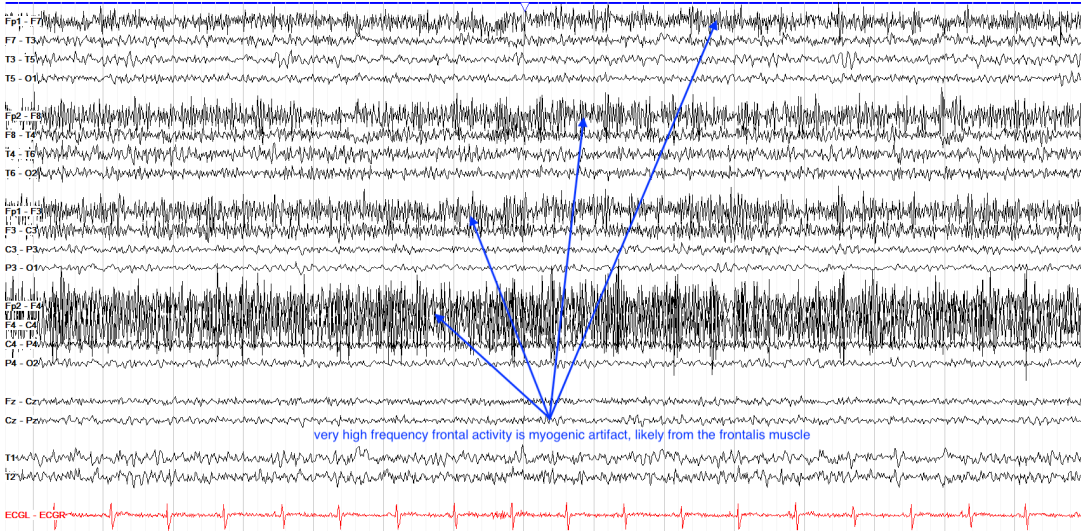


Figure 7.1: Example of artifacts in EEG data detected via visual inspection.

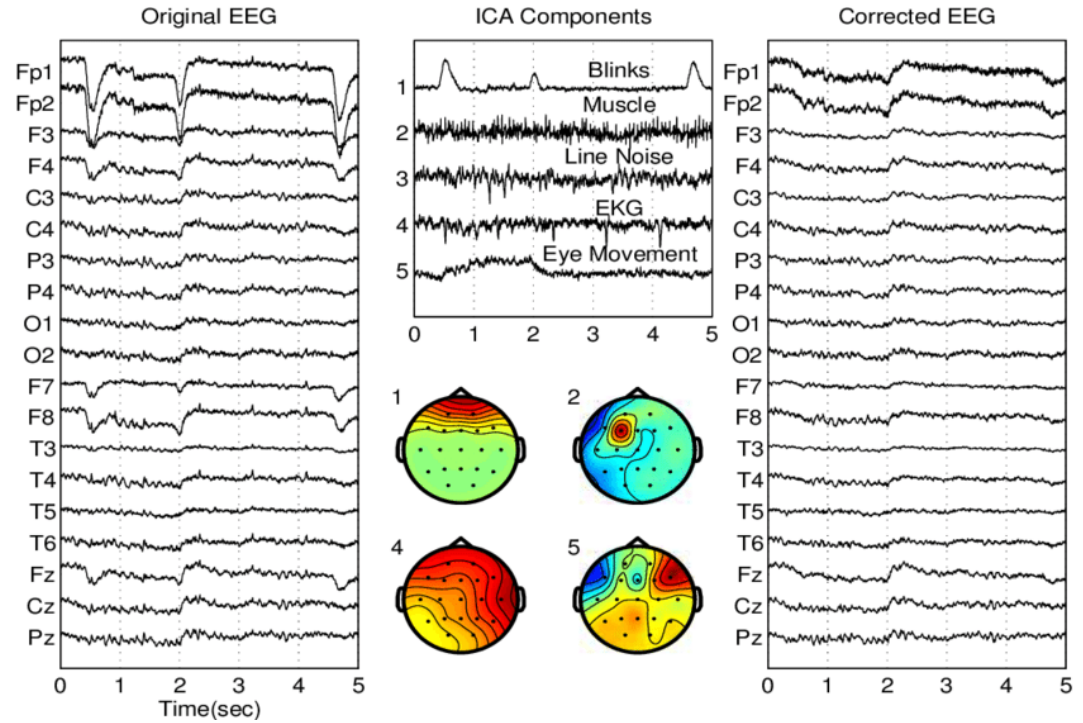


Figure 7.2: Visualization of ICA components.

### 7.3.1 Quality Metrics

Quality metrics for evaluating preprocessing include:

- **Signal-to-Noise Ratio (SNR):** Measures the relative strength of the neural signal to background noise.
- **Power Spectral Density (PSD):** Evaluates changes in frequency-domain characteristics before and after artifact removal.
- **Artifact Detection Success Rate:** Quantifies the proportion of artifacts correctly identified and removed.

### 7.3.2 Comparing Pre- and Post-Correction Data

Visual and statistical comparisons of pre- and post-correction data can highlight the effectiveness of artifact removal. In MNE-Python, data before and after correction can be plotted side by side:

```
# Plot pre- and post-correction data  
raw.plot(title="Before_Correction")  
raw_ica.plot(title="After_ICA_Correction")
```

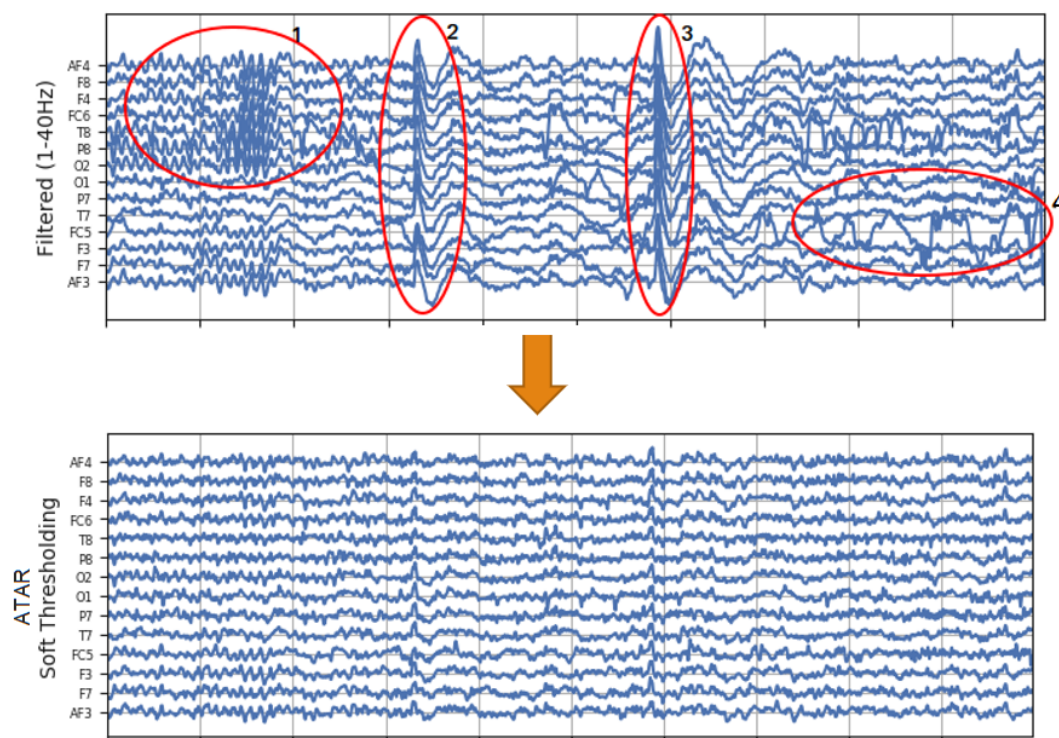


Figure 7.3: Comparison of EEG data before and after artifact correction.

# Bibliography

- [1] Gramfort, A., et al. (2013). MEG and EEG data analysis with MNE-Python. *Frontiers in Neuroscience*, 7, 267.
- [2] MNE-Python Documentation. (2023). Retrieved from <https://mne.tools/stable/index.html>.
- [3] Cohen, M. X. (2014). *Analyzing Neural Time Series Data: Theory and Practice*. MIT Press.
- [4] Delorme, A., Makeig, S. (2004). EEGLAB: An open-source toolbox for analysis of single-trial EEG dynamics. *Journal of Neuroscience Methods*, 134(1), 9–21.

## **Part V**

# **Advanced EEG Signal Analysis**

# Chapter 8

## Feature Extraction and Time-Frequency Analysis

Analyzing EEG data involves extracting meaningful features that characterize brain activity. These features can be derived from different domains, including time, frequency, and time-frequency representations. This chapter explores key methods for feature extraction and analysis in EEG data.

### 8.1 Time-Domain Analysis

Time-domain analysis focuses on features derived directly from the amplitude of EEG signals over time.

#### 8.1.1 Event-Related Potentials (ERPs)

Event-Related Potentials (ERPs) are time-locked neural responses to specific stimuli or events. They are often averaged across multiple trials to enhance signal-to-noise ratio.

- **Components:** Common ERP components include the P300 (positive deflection around 300 ms) and N170 (negative deflection associated with facial recognition).
- **Applications:** Used in cognitive neuroscience to study attention, perception, and decision-making.

In MNE-Python, ERPs are computed as follows:



```
# Create epochs
epochs = mne.Epochs(raw, events, event_id={'stimulus': 1}, tmin=-0.2, tmax=0.5, b

# Compute evoked response
evoked = epochs.average()

# Plot ERP
evoked.plot()
```

### 8.1.2 Microstate Analysis

Microstates are brief, stable topographical patterns in EEG signals that reflect global brain states.

- **Duration:** Typically last 50–120 ms.
- **Applications:** Associated with cognitive and clinical states, including schizophrenia and Alzheimer's disease.

MNE-Python supports microstate analysis through additional packages or custom implementations.

## 8.2 Frequency-Domain Analysis

Frequency-domain analysis involves studying the power and phase of EEG signals in specific frequency bands, such as delta, theta, alpha, beta, and gamma.

### 8.2.1 Power Spectral Density (PSD)

The Power Spectral Density (PSD) represents the distribution of power across frequencies and provides insights into the dominant rhythms of brain activity.

```
# Compute and plot PSD
raw.plot_psd(fmin=1, fmax=40)
```

### 8.2.2 Fast Fourier Transform (FFT) Techniques

The Fast Fourier Transform (FFT) is a computationally efficient algorithm for converting time-domain signals into the frequency domain.

```
# Compute FFT
import numpy as np
from scipy.fftpack import fft

data = raw.get_data() # Extract EEG data
fft_data = fft(data, axis=1) # Perform FFT along time axis
```

FFT is useful for identifying the dominant frequencies in EEG signals, aiding in band-specific analyses.

## 8.3 Time-Frequency Representations

Time-frequency analysis combines temporal and spectral information, enabling the study of how frequency components change over time.

### 8.3.1 Short-Time Fourier Transform (STFT)

The Short-Time Fourier Transform (STFT) divides the signal into overlapping time windows and computes the FFT for each window.

```
# Perform STFT using MNE-Python
from mne.time_frequency import tfr_multitaper

# Compute time-frequency representation
tfr = tfr_multitaper(epochs, fmin=1, fmax=40, n_cycles=5, time_bandwidth=2.0)
tfr.plot()
```

### 8.3.2 Wavelet Transformations

Wavelet transformations use wavelets, localized functions in time and frequency, to decompose signals.

```
# Perform wavelet transformation
from mne.time_frequency import tfr_morlet

# Compute wavelet transform
frequencies = np.arange(1, 40, 1) # Define frequency range
tfr_wavelet = tfr_morlet(epochs, freqs=frequencies, n_cycles=frequencies / 2)
tfr_wavelet.plot()
```

Wavelets are especially effective for analyzing non-stationary signals like EEG.

## 8.4 Spatial Analysis

Spatial analysis involves examining the topographical distribution of EEG activity across the scalp.

### 8.4.1 Source Localization

Source localization estimates the cortical sources generating the observed EEG signals.

- **Methods:** Includes dipole fitting, beamforming, and minimum norm estimation.
- **Applications:** Used in epilepsy research and brain-computer interfaces.

```
# Source localization example (requires forward model and inverse operator)
from mne.minimum_norm import apply_inverse

stc = apply_inverse(evoked, inverse_operator, lambda2=1.0 / 9.0, method='dSPM')
stc.plot()
```

### 8.4.2 Topographical Mapping

Topographical maps visualize the spatial distribution of EEG activity across the scalp.

```
# Plot topographical map
evoked.plot_topomap(times=[0.1, 0.2, 0.3], ch_type='eeg')
```

These advanced techniques provide comprehensive insights into the temporal, spectral, and spatial characteristics of EEG data, supporting diverse research and clinical applications.

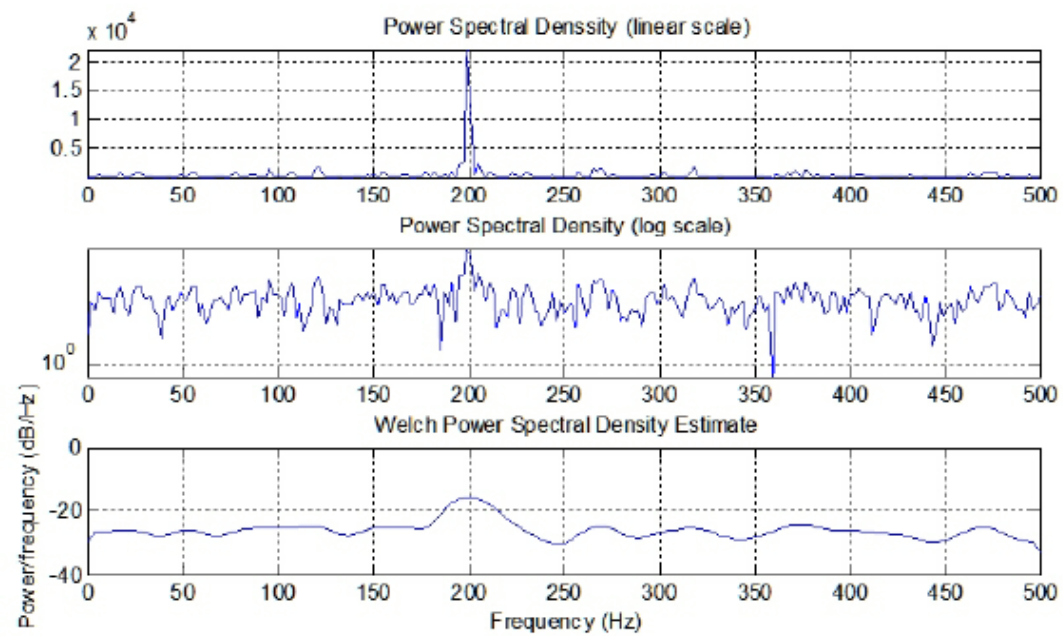


Figure 8.1: Power Spectral Density of EEG data, highlighting dominant frequency bands.

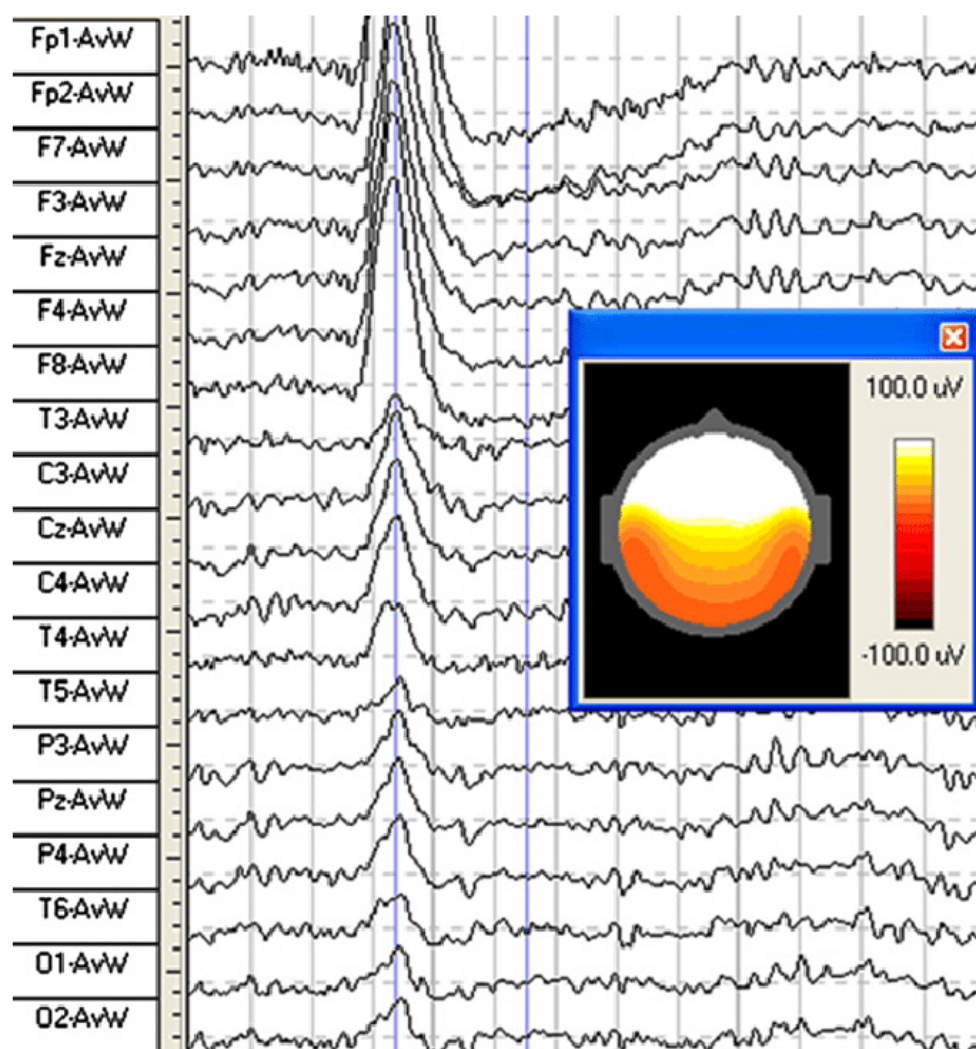


Figure 8.2: Topographical map of EEG activity at specific time points.

# Bibliography

- [1] Gramfort, A., et al. (2013). MEG and EEG data analysis with MNE-Python. *Frontiers in Neuroscience*, 7, 267.
- [2] MNE-Python Documentation. (2023). Retrieved from <https://mne.tools/stable/index.html>.
- [3] Cohen, M. X. (2014). *Analyzing Neural Time Series Data: Theory and Practice*. MIT Press.
- [4] Delorme, A., Makeig, S. (2004). EEGLAB: An open-source toolbox for analysis of single-trial EEG dynamics. *Journal of Neuroscience Methods*, 134(1), 9–21.

# Chapter 9

## Machine Learning and EEG

The integration of machine learning (ML) techniques with electroencephalography (EEG) data analysis has significantly advanced the understanding and interpretation of neural signals. This chapter provides an overview of fundamental ML concepts, feature selection methods, classification algorithms, deep learning approaches, and practical implementation strategies using Python.

### 9.1 Introduction to Machine Learning Concepts

Machine learning encompasses algorithms that enable computers to learn patterns from data and make decisions or predictions without explicit programming.

#### 9.1.1 Supervised vs. Unsupervised Learning

- **Supervised Learning:** Involves training models on labeled datasets, where the input data is paired with the correct output. Common applications include classification and regression tasks. For instance, classifying EEG signals into different mental states based on labeled examples.
- **Unsupervised Learning:** Deals with unlabeled data, aiming to uncover hidden structures or patterns. Techniques such as clustering and dimensionality reduction fall under this category. An example is identifying distinct patterns in EEG data without predefined labels.

### 9.1.2 Model Evaluation Metrics

Evaluating the performance of ML models is crucial. Common metrics include:

- **Accuracy:** The proportion of correctly predicted instances among the total instances.
- **Precision:** The ratio of true positive predictions to the total predicted positives.
- **Recall (Sensitivity):** The ratio of true positive predictions to all actual positives.
- **F1-Score:** The harmonic mean of precision and recall, providing a balance between the two.
- **Area Under the Receiver Operating Characteristic Curve (AUC-ROC):** Measures the ability of the model to distinguish between classes.

These metrics provide insights into the model's performance, guiding improvements and ensuring reliability.

## 9.2 Feature Selection and Dimensionality Reduction

EEG data often comprises high-dimensional features, necessitating techniques to reduce dimensionality while preserving essential information.

### 9.2.1 Principal Component Analysis (PCA)

PCA is a linear dimensionality reduction technique that transforms data into a set of orthogonal components, capturing the maximum variance.

- **Application:** Reduces the number of features in EEG data, mitigating the curse of dimensionality and enhancing computational efficiency.
- **Implementation:**

```
from sklearn.decomposition import PCA

pca = PCA(n_components=10)
reduced_data = pca.fit_transform(eeg_data)
```



### 9.2.2 t-Distributed Stochastic Neighbor Embedding (t-SNE)

t-SNE is a non-linear dimensionality reduction technique that visualizes high-dimensional data by mapping it into a low-dimensional space.

- **Application:** Visualizes complex patterns in EEG data, aiding in the identification of clusters or anomalies.
- **Implementation:**

```
from sklearn.manifold import TSNE

tsne = TSNE(n_components=2)
embedded_data = tsne.fit_transform(eeg_data)
```

## 9.3 Classification Algorithms

Various ML algorithms are employed to classify EEG signals into distinct categories.

### 9.3.1 Support Vector Machines (SVM)

SVMs are supervised learning models that find the optimal hyperplane separating classes in the feature space.

- **Application:** Effective in classifying EEG signals, especially when the data is linearly separable.
- **Implementation:**

```
from sklearn.svm import SVC

svm = SVC(kernel='linear')
svm.fit(X_train, y_train)
predictions = svm.predict(X_test)
```

### 9.3.2 Random Forests

Random Forests are ensemble learning methods that construct multiple decision trees and aggregate their outputs.

- **Application:** Handle high-dimensional EEG data and provide insights into feature importance.
- **Implementation:**

```
from sklearn.ensemble import RandomForestClassifier

rf = RandomForestClassifier(n_estimators=100)
rf.fit(X_train, y_train)
predictions = rf.predict(X_test)
```

### 9.3.3 K-Nearest Neighbors (KNN)

KNN is a simple, instance-based learning algorithm that classifies data points based on the majority label among their k-nearest neighbors.

- **Application:** Useful for EEG classification tasks with well-defined clusters.
- **Implementation:**

```
from sklearn.neighbors import KNeighborsClassifier

knn = KNeighborsClassifier(n_neighbors=5)
knn.fit(X_train, y_train)
predictions = knn.predict(X_test)
```

## 9.4 Deep Learning Approaches

Deep learning models, particularly neural networks, have shown promise in capturing complex patterns in EEG data.

### 9.4.1 Convolutional Neural Networks (CNN) for EEG

CNNs are designed to process grid-like data structures and are effective in capturing spatial hierarchies.

- **Application:** Applied to EEG data by treating the signal as a 2D image, capturing spatial and temporal features.

- **Implementation:**

```

from tensorflow.keras.models import Sequential
from tensorflow.keras.layers import Conv2D, MaxPooling2D, Flatten, Dense

model = Sequential([
    Conv2D(32, kernel_size=(3, 3), activation='relu', input_shape=input_shape),
    MaxPooling2D(pool_size=(2, 2)),
    Flatten(),
    Dense(128, activation='relu'),
    Dense(num_classes, activation='softmax')
])
model.compile(optimizer='adam', loss='categorical_crossentropy', metrics=
model.fit(X_train, y_train, epochs=10, batch_size=32)

```

### 9.4.2 Recurrent Neural Networks (RNN) and LSTM

Recurrent Neural Networks (RNNs) are designed to process sequential data by maintaining a memory of previous inputs, making them suitable for time-series analysis such as EEG signals. However, standard RNNs can suffer from the vanishing gradient problem, hindering their ability to learn long-term dependencies. Long Short-Term Memory (LSTM) units address this issue by incorporating mechanisms to retain information over extended periods.

- **Application:** LSTMs are effective in modeling temporal dependencies in EEG data, enabling the analysis of dynamic brain activities and aiding in tasks like seizure detection and mental state classification.

- **Implementation:**

```

from tensorflow.keras.models import Sequential
from tensorflow.keras.layers import LSTM, Dense

# Define the model
model = Sequential([
    LSTM(100, activation='relu', input_shape=(n_timesteps, n_features)),
    Dense(1)
])

```

```

# Compile the model
model.compile(optimizer='adam', loss='mse')

# Fit the model
model.fit(X_train, y_train, epochs=10, batch_size=32)

```

In this example, ' $n_{timesteps}$ ' represents the number of timesteps in each EEG sequence, and ' $n_{features}$ ' denotes the

## 9.5 Model Implementation with Python

Implementing machine learning models for EEG analysis in Python involves utilizing libraries tailored for both traditional machine learning and deep learning approaches.

### 9.5.1 Scikit-Learn for Traditional ML

Scikit-Learn is a versatile library that provides simple and efficient tools for data analysis and modeling.

- **Application:** Suitable for implementing traditional machine learning algorithms such as Support Vector Machines (SVM), Random Forests, and K-Nearest Neighbors (KNN) for EEG data classification and regression tasks.
- **Implementation:**

```

from sklearn.ensemble import RandomForestClassifier
from sklearn.metrics import accuracy_score

# Initialize the model
rf = RandomForestClassifier(n_estimators=100)

# Train the model
rf.fit(X_train, y_train)

# Make predictions
predictions = rf.predict(X_test)

# Evaluate the model
accuracy = accuracy_score(y_test, predictions)

```

```
print(f 'Accuracy: {accuracy:.2f}')
```

In this example, a Random Forest classifier is trained on EEG data to predict class labels. The model's performance is evaluated using accuracy, but other metrics like precision, recall, and F1-score can also be employed depending on the specific application.

### 9.5.2 TensorFlow and PyTorch for Deep Learning

TensorFlow and PyTorch are prominent deep learning frameworks that facilitate the development and training of complex neural networks.

- **TensorFlow:**

- **Application:** TensorFlow, often used with its high-level Keras API, simplifies the construction of deep learning models, including CNNs and RNNs, for EEG data analysis.
- **Implementation:**

```
from tensorflow.keras.models import Sequential
from tensorflow.keras.layers import Conv2D, MaxPooling2D, Flatten,

# Define the model
model = Sequential([
    Conv2D(32, kernel_size=(3, 3), activation='relu', input_shape=(
    MaxPooling2D(pool_size=(2, 2)),
    Flatten(),
    Dense(128, activation='relu'),
    Dense(num_classes, activation='softmax')
])

# Compile the model
model.compile(optimizer='adam', loss='categorical_crossentropy', m

# Train the model
model.fit(X_train, y_train, epochs=10, batch_size=32)
```

This example demonstrates a Convolutional Neural Network (CNN) implemented using TensorFlow's Keras API, designed for classifying EEG data into multiple categories.

- **PyTorch:**

- **Application:** PyTorch offers dynamic computation graphs and an intuitive interface, making it suitable for research and applications requiring flexibility in model design, such as custom neural network architectures for EEG analysis.

- **Implementation:**

```
import torch
import torch.nn as nn
import torch.optim as optim
from torch.utils.data import DataLoader, TensorDataset

# Define the model
class EEGNet(nn.Module):
    def __init__(self):
        super(EEGNet, self).__init__()
        self.conv1 = nn.Conv2d(1, 32, kernel_size=(3, 3))
        self.pool = nn.MaxPool2d(kernel_size=(2, 2))
        self.fc1 = nn.Linear(32 * 6 * 6, 128)
        self.fc2 = nn.Linear(128, num_classes)

    def forward(self, x):
        x = self.pool(F.relu(self.conv1(x)))
        x = x.view(-1, 32 * 6 * 6)
        x = F.relu(self.fc1(x))
        x = self.fc2(x)
        return x

# Initialize the model, loss function, and optimizer
model = EEGNet()
criterion = nn.CrossEntropyLoss()
optimizer = optim.Adam(model.parameters(), lr=0.001)

# Prepare data loaders
train_dataset = TensorDataset(torch.tensor(X_train).float(), torch.tensor(y_train))
train_loader = DataLoader(train_dataset, batch_size=32, shuffle=True)

# Train the model
```

```
for epoch in range(10):  
    for inputs, labels in train_loader:  
        optimizer.zero_grad()  
        outputs = model(inputs)  
        loss = criterion(outputs, labels)  
        loss.backward()  
        optimizer.step()
```

This example illustrates a simple CNN implemented in PyTorch for EEG data classification. The model is trained using a DataLoader to handle batches of data efficiently.

Both TensorFlow and PyTorch provide robust tools for developing deep learning models tailored to the complexities of EEG data, enabling advancements in brain-computer interfaces, cognitive state monitoring, and neurological disorder detection.

# Chapter 10

## Visualizing EEG Data

Effective visualization of EEG data is crucial for interpreting complex neural signals and identifying patterns. This chapter explores various techniques and tools for visualizing EEG data, including time-series plots, spectral and topographical representations, advanced methods, and custom visualizations using Python libraries.

### 10.1 Time-Series Visualization

Time-series plots display EEG signals over time, providing insights into temporal patterns and anomalies.

#### 10.1.1 Multi-Channel Plotting

EEG recordings typically involve multiple channels corresponding to different scalp locations. Visualizing all channels simultaneously aids in comparing activity across regions.

- **Application:** Detecting artifacts, identifying event-related potentials, and assessing overall signal quality.
- **Implementation:**

```
import mne
```

```
# Load raw EEG data
```

```
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)
```



```
# Plot all channels  
raw.plot(n_channels=64, scalings='auto', title='EEG_Data', show=True)
```

This code utilizes MNE-Python to load and plot EEG data, displaying all channels with automatic scaling for optimal visualization.

### 10.1.2 Interactive Visualization Tools

Interactive tools enhance data exploration by allowing zooming, panning, and channel selection.

- **Application:** Facilitating detailed examination of specific time windows or channels.
- **Implementation:**

```
import mne  
  
# Load raw EEG data  
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)  
  
# Launch interactive plot  
raw.plot(block=True)
```

Setting 'block=True' enables interactive features, allowing users to navigate through the data effectively.

## 10.2 Spectral and Topographical Plots

These plots represent the frequency content and spatial distribution of EEG signals.

### 10.2.1 Spectrograms and Heatmaps

Spectrograms display how the spectral content of a signal evolves over time, while heatmaps represent data values across two dimensions.

- **Application:** Identifying frequency-specific patterns, such as alpha or beta rhythms, and their temporal dynamics.
- **Implementation:**

```

import mne
import matplotlib.pyplot as plt

# Load raw EEG data
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)

# Plot spectrogram for a specific channel
raw.plot_psd(picks=['EEG_001'], fmin=0, fmax=50, show=True)

```

This code generates a power spectral density (PSD) plot for a selected EEG channel, highlighting its frequency components.

### 10.2.2 Scalp Maps and 3D Brain Representations

Scalp maps visualize the spatial distribution of EEG activity across the scalp, and 3D representations provide insights into cortical sources.

- **Application:** Localizing neural activity and understanding spatial patterns.
- **Implementation:**

```

import mne

# Load evoked data
evoked = mne.read_evokeds('sample-ave.fif', condition='Auditory', baseline

# Plot topomap at a specific time point
evoked.plot_topomap(times=0.1, ch_type='eeg', show=True)

```

This script plots a topographical map of EEG activity at 100 ms post-stimulus, illustrating the spatial distribution of the response.

## 10.3 Advanced Visualization Techniques

Advanced methods provide deeper insights into neural connectivity and dynamics.

### 10.3.1 Connectivity Graphs

Connectivity graphs depict functional connections between different brain regions based on EEG data.

- **Application:** Studying neural networks and communication pathways.
- **Implementation:**

```
import mne
from mne.connectivity import spectral_connectivity

# Load raw EEG data
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)

# Compute connectivity
con, freqs, times, n_epochs, n_tapers = spectral_connectivity(
    raw, method='pli', sfreq=raw.info['sfreq'], fmin=8, fmax=13, faverage=

# Plot connectivity circle
mne.viz.plot_connectivity_circle(con[:, :, 0], raw.info['ch_names'], title=
```

This code computes phase lag index (PLI) connectivity in the alpha band and visualizes it as a circular graph, highlighting connections between channels.

### 10.3.2 Dynamic Visualizations and Animations

Animating EEG data over time can reveal dynamic changes in brain activity.

- **Application:** Observing temporal evolution of neural responses and identifying transient events.
- **Implementation:**

```
import mne

# Load raw EEG data
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)

# Create an animation of the EEG data
raw.animate(duration=10, blit=True)
```

The ‘animate’ function in MNE-Python generates an animation of the EEG signals over a specified duration, facilitating the observation of dynamic patterns.

## 10.4 Custom Visualization with Matplotlib and MNE-Python

For tailored visualizations, combining Matplotlib with MNE-Python offers flexibility.

- **Application:** Creating customized plots to highlight specific aspects of EEG data.
- **Implementation:**

```
import mne
import matplotlib.pyplot as plt

# Load raw EEG data
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)

# Select a specific channel
data, times = raw['EEG_001']

# Plot the data
plt.figure()
plt.plot(times, data.T)
plt.title('EEG_Channel_001')
plt.xlabel('Time(s)')
plt.ylabel('Amplitude( V )')
plt.show()
```

This script extracts data from a specific EEG channel and plots it using Matplotlib, allowing for customized styling and annotations.

# Appendix A

## Installation Guides and Troubleshooting

### A.1 Setting Up Python and MNE-Python

To effectively analyze EEG data, it's essential to set up a Python environment with the MNE-Python library. Follow these steps:

1. **Install Python:** Ensure that Python 3.9 or higher is installed on your system. You can download it from the official Python website.
2. **Install MNE-Python:** Use pip to install MNE-Python with minimal dependencies:

```
pip install --upgrade mne
```

For more detailed instructions, refer to the MNE-Python installation guide. :contentReference[oaicite:0]index=0

3. **Verify Installation:** After installation, verify by importing MNE in a Python script or interactive session:

```
import mne
print(mne.__version__)
```

This should display the installed MNE-Python version.

### A.2 Common Issues and Solutions

While setting up and using MNE-Python, you might encounter some common issues:

- **Issue:** ImportError: No module named 'mne'
  - **Solution:** Ensure that MNE-Python is installed in the correct Python environment. Activate the environment where MNE is installed or install MNE in the active environment.
- **Issue:** Compatibility problems with other packages
  - **Solution:** Use virtual environments to manage dependencies and avoid conflicts. Tools like venv or conda can help create isolated environments.
- **Issue:** Visualization errors or missing GUI support
  - **Solution:** Ensure that the necessary visualization backends are installed. For example, installing PyQt5 can resolve GUI-related issues:

```
pip install pyqt5
```

For more troubleshooting tips, consult the MNE-Python documentation.

# Appendix B

## Sample Datasets and Code Examples

### B.1 Accessing Public EEG Databases

Utilizing public EEG datasets is invaluable for practice and research. Notable repositories include:

- **OpenNeuro:** An open platform for sharing neuroimaging data, including EEG datasets.
- **EEG-Datasets on GitHub:** A curated list of public EEG datasets suitable for various applications.
- **PhysioNet:** Offers a vast collection of physiological signal databases, including EEG data.

### B.2 Annotated Code Snippets for Key Procedures

Here are some essential code snippets for EEG data analysis using MNE-Python:

#### B.2.1 Loading and Plotting Raw EEG Data

```
import mne
```

```
# Load raw EEG data
```

```
raw = mne.io.read_raw_fif('sample_data.fif', preload=True)
```

```
# Plot raw data
```

```
raw.plot(n_channels=10, duration=5, scalings='auto')
```

### B.2.2 Filtering EEG Data

```
# Apply a band-pass filter from 1 to 40 Hz  
raw.filter(l_freq=1.0, h_freq=40.0)
```

### B.2.3 Epoching Data Around Events

Epoching refers to segmenting continuous EEG data into time windows around specific events, such as stimuli or responses. These segments, called epochs, facilitate analysis of time-locked neural activity.

```
# Define events  
events = mne.find_events(raw)  
  
# Create epochs around events with specific event ID  
epochs = mne.Epochs(raw, events, event_id={'stimulus': 1},  
                    tmin=-0.2, tmax=0.5, baseline=(None, 0))  
  
# Plot epoch data  
epochs.plot()
```

The ‘tmin’ and ‘tmax’ parameters define the time range relative to the event (e.g., -200 ms to 500 ms), while the ‘baseline’ argument corrects for baseline drift.

### B.2.4 Computing and Plotting Evoked Responses

```
# Compute evoked response  
evoked = epochs.average()  
  
# Plot evoked response  
evoked.plot()
```



# Appendix C

## Mathematical Foundations

### C.1 Linear Algebra and Signal Processing Basics

Understanding the following concepts is crucial for EEG data analysis:

- **Fourier Transform:** Decomposes signals into their constituent frequencies, aiding in frequency domain analysis.
- **Convolution:** Describes how the shape of one function is modified by another, fundamental in filtering operations.
- **Eigenvalues and Eigenvectors:** Key in dimensionality reduction techniques like Principal Component Analysis (PCA).

### C.2 Statistical Methods in EEG Analysis

Statistical techniques are essential for interpreting EEG data:

- **Hypothesis Testing:** Determines the significance of observed effects, such as differences between conditions.
- **Correlation Analysis:** Assesses relationships between EEG signals and behavioral measures.
- **Time-Frequency Analysis:** Examines how spectral content evolves over time, revealing dynamic neural processes.

# Appendix D

## Additional Resources

### D.1 Recommended Books and Articles

- **"An Introduction to the Event-Related Potential Technique"** by Steven J. Luck: A comprehensive guide to ERP methods and applications.
- **"Handbook of EEG Interpretation"** by William O. Tatum: Offers insights into EEG patterns and their clinical significance.
- **"EEG Signal Processing"** by Saeid Sanei and J.A. Chambers: Explores advanced signal processing techniques for EEG analysis.

### D.2 Online Courses and Tutorials

- **Coursera's "EEG Signal Processing"**: Provides foundational knowledge in EEG analysis techniques.
- **MNE-Python Tutorials**: Offers practical guides and examples for EEG data analysis using MNE-Python.

# Appendix E

## Glossary of Terms

- **EEG (Electroencephalography):** A method to record electrical activity of the brain via electrodes placed on the scalp.
- **Epoch:** A segment of EEG data time-locked to a specific event, used in event-related analyses.
- **Artifact:** Unwanted interference or noise in EEG data, often from sources like muscle movements

# Appendix F

## Practice Project: Generating Synthetic EEG Data and Implementing Clustering Algorithms

This project walks through the steps of generating synthetic EEG data, preprocessing it using a bandpass filter, clustering the data using K-means, and visualizing the results. The goal is to provide a hands-on approach to understanding EEG data manipulation and machine learning techniques.

### F.1 Step 1: Importing Necessary Libraries

To begin, import the necessary libraries for numerical operations, signal processing, machine learning, and visualization:

```
import numpy as np
from scipy.signal import butter, lfilter
from sklearn.cluster import KMeans
from sklearn.decomposition import PCA
import matplotlib.pyplot as plt
```

#### Explanation:

- `numpy`: Provides support for large, multi-dimensional arrays and matrices.
- `scipy.signal`: Contains functions for filtering and signal processing.
- `sklearn.cluster` and `sklearn.decomposition`: Implements machine learning algorithms for clustering and dimensionality reduction.

- `matplotlib.pyplot`: Enables visualization of data and clustering results.

## F.2 Step 2: Generating Synthetic EEG Data

Define a function to generate synthetic EEG data. Each "electrode" simulates a sine wave with added noise:

```
def generate_eeg_data(samples=1000, electrodes=10):
    freq = np.random.uniform(0.5, 30, electrodes)
    # Random frequencies
    time = np.linspace(0, 1, samples)
    eeg_data = np.array([np.sin(2 * np.pi * f * time) + 0.5 * np.random.randn()
                        for f in freq])
    return eeg_data, time
```

### Explanation:

- `freq`: Assigns random frequencies between 0.5 and 30 Hz to each electrode, simulating different neural rhythms.
- `time`: Creates a time vector for one second of data.
- `eeg_data`: Generates sine waves with added Gaussian noise, simulating noisy neural signals.

**Output:** A matrix where rows represent electrode signals and columns represent time points.

## F.3 Step 3: Applying a Bandpass Filter

Define a function to filter the EEG data. A bandpass filter allows only signals within a specified frequency range to pass through:

```
def bandpass_filter(data, lowcut, highcut, fs, order=5):
    nyquist = 0.5 * fs
    low = lowcut / nyquist
    high = highcut / nyquist
    b, a = butter(order, [low, high], btype='band')
    filtered_data = lfilter(b, a, data, axis=1)
    return filtered_data
```

### Explanation:

- `lowcut` and `highcut`: Define the passband range in Hz.
- `fs`: Sampling frequency.
- `butter`: Designs a Butterworth filter, which is smooth in the frequency domain.
- `lfilter`: Applies the filter to the EEG data along the time dimension (`axis=1`).

## F.4 Step 4: Clustering the EEG Data

Apply K-means clustering to the filtered data to group the signals into clusters based on their patterns:

```
from sklearn.cluster import KMeans

reshaped_data = filtered_data.reshape(filtered_data.shape[0], -1)
kmeans = KMeans(n_clusters=n_clusters, random_state=42)
clusters = kmeans.fit_predict(reshaped_data)
```

### Explanation:

- `reshape`: Converts the data into a 2D array where each row represents an electrode's signal over time.
- `KMeans`: Initializes the clustering algorithm with a specified number of clusters (`n_clusters`).
- `fit_predict`: Fits the clustering model to the data and assigns cluster labels to each signal.

## F.5 Step 5: Visualizing the Results

Use PCA to reduce the dimensionality of the data and visualize the clusters in a 2D space. Also, plot the raw and filtered signals:

```
from sklearn.decomposition import PCA

pca = PCA(n_components=2)
pca_results = pca.fit_transform(reshaped_data)

# Visualization code (abridged)
plt.scatter(pca_results[:, 0], pca_results[:, 1], c=clusters)
```

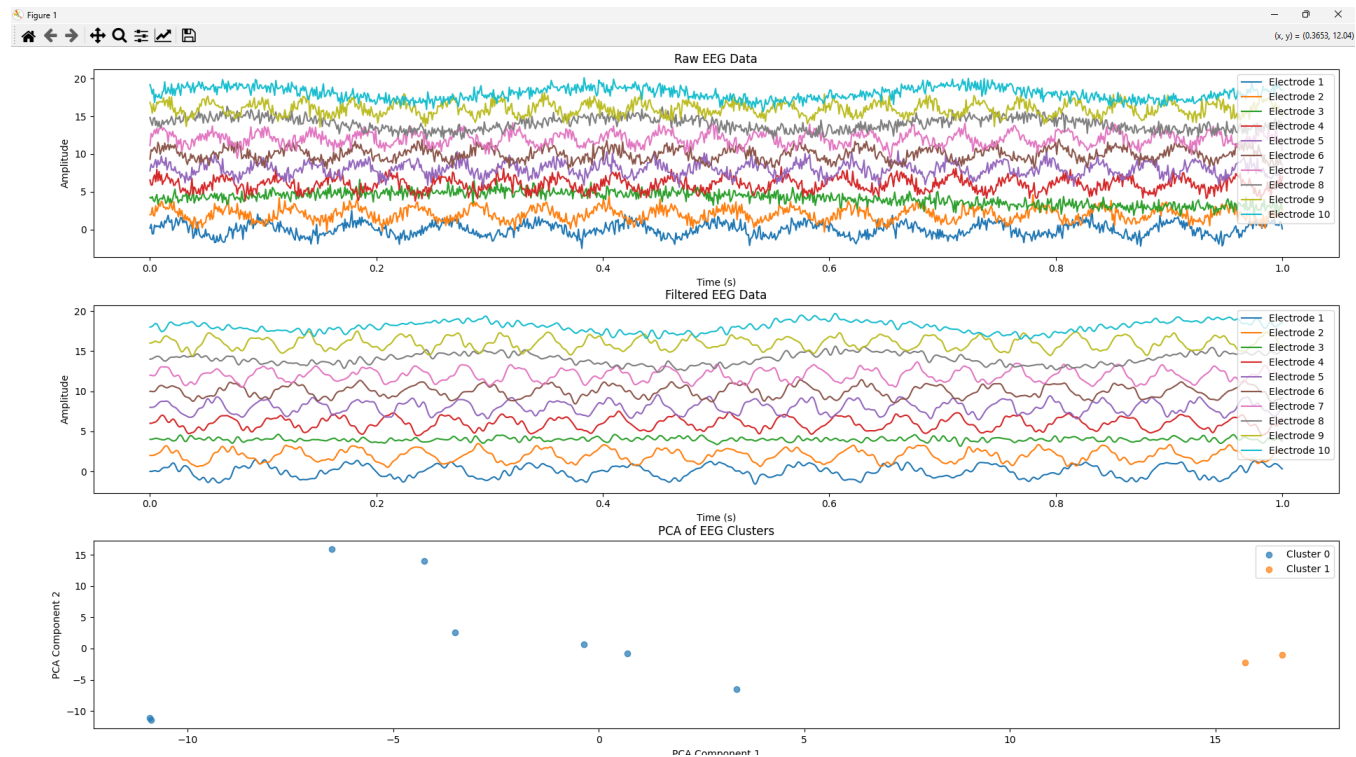


Figure F.1: Visualized EEG signal, pre and post-filtering, and 2D PCA

```
plt.xlabel('PCA_Component_1')
plt.ylabel('PCA_Component_2')
plt.title('PCA of EEG Clusters')
plt.show()
```


**Explanation:**

- **PCA**: Reduces the data to two principal components, simplifying visualization.
- **scatter**: Plots the clustered data in the reduced space.
- **c=clusters**: Colors the data points based on their cluster assignments.

## F.6 Step 6: Running the Program

Finally, combine the functions in a main program that prompts the user for parameters and runs the entire pipeline:

```
if __name__ == "__main__":
    samples = int(input("Enter the number of samples: "))
    lowcut = float(input("Enter the low cutoff frequency: "))
```

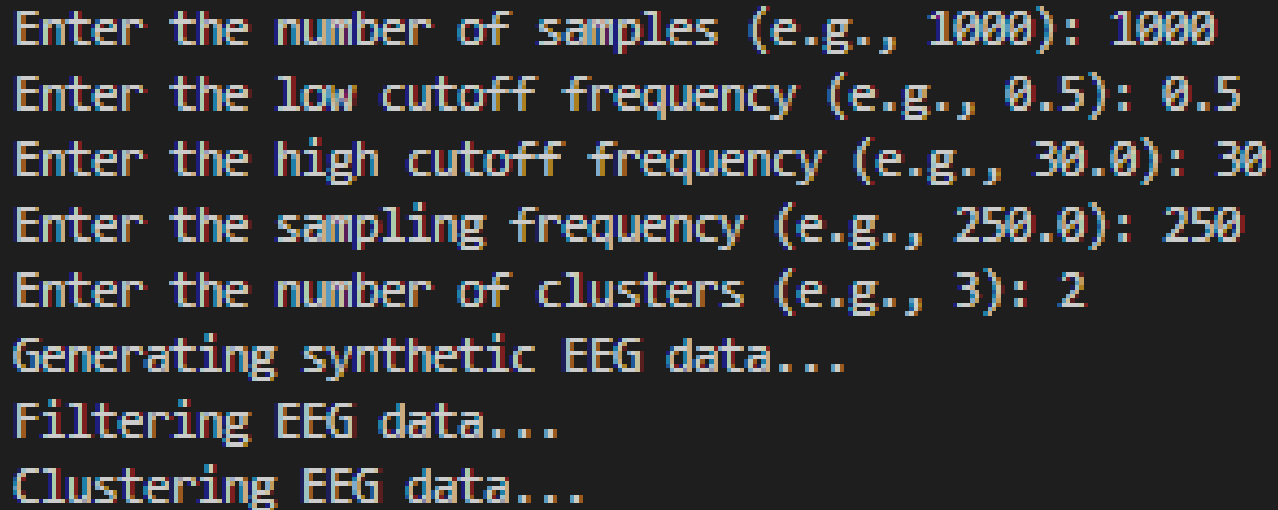A terminal window with a black background and light blue text. The text shows a series of prompts and user inputs: 'Enter the number of samples (e.g., 1000): 1000', 'Enter the low cutoff frequency (e.g., 0.5): 0.5', 'Enter the high cutoff frequency (e.g., 30.0): 30', 'Enter the sampling frequency (e.g., 250.0): 250', 'Enter the number of clusters (e.g., 3): 2', followed by three status messages: 'Generating synthetic EEG data...', 'Filtering EEG data...', and 'Clustering EEG data...'.

```
Enter the number of samples (e.g., 1000): 1000
Enter the low cutoff frequency (e.g., 0.5): 0.5
Enter the high cutoff frequency (e.g., 30.0): 30
Enter the sampling frequency (e.g., 250.0): 250
Enter the number of clusters (e.g., 3): 2
Generating synthetic EEG data...
Filtering EEG data...
Clustering EEG data...
```

Figure F.2: Running the program within the terminal

```
highcut = float(input("Enter the high cutoff frequency: "))
fs = float(input("Enter the sampling frequency: "))
n_clusters = int(input("Enter the number of clusters: "))

eeg_data, time = generate_eeg_data(samples=samples)
filtered_data = bandpass_filter(eeg_data, lowcut, highcut, fs)
# Additional processing and visualization here...
```

This section ties together data generation, filtering, clustering, and visualization into a cohesive workflow.



# Bibliography

- [1] Gramfort, A., Luessi, M., Larson, E., Engemann, D. A., Strohmeier, D., Brodbeck, C., ... & Hämäläinen, M. S. (2013). MNE software for processing MEG and EEG data. *NeuroImage*, 86, 446-460.
- [2] Luck, S. J. (2014). *An Introduction to the Event-Related Potential Technique* (2nd ed.). MIT Press.
- [3] Tatum, W. O. (2014). *Handbook of EEG Interpretation* (2nd ed.). Demos Medical Publishing.
- [4] Sanei, S., & Chambers, J. A. (2007). *EEG Signal Processing*. Wiley-Interscience.
- [5] Poldrack, R. A., Gorgolewski, K. J., & Esteban, O. (2017). OpenNeuro: A free online platform for sharing and analysis of neuroimaging data. *Organization for Human Brain Mapping*.
- [6] Roy, Y., Banville, H., Albuquerque, I., Gramfort, A., Falk, T. H., & Faubert, J. (2019). Deep learning-based electroencephalography analysis: A systematic review. *Journal of Neural Engineering*, 16(5), 051001.
- [7] Goldberger, A. L., Amaral, L. A. N., Glass, L., Hausdorff, J. M., Ivanov, P. C., Mark, R. G., ... & Stanley, H. E. (2000). PhysioBank, PhysioToolkit, and PhysioNet: Components of a new research resource for complex physiologic signals. *Circulation*, 101(23), e215-e220.
- [8] Makeig, S., & Delorme, A. (2012). EEG Signal Processing. *Coursera*.
- [9] Gramfort, A., Luessi, M., Larson, E., Engemann, D. A., Strohmeier, D., Brodbeck, C., ... & Hämäläinen, M. S. (2014). MNE-Python: A comprehensive open-source environment for EEG and MEG data analysis. *Frontiers in Neuroscience*, 7, 267.
- [10] Acharya, U. R., Sree, S. V., Swapna, G., Martis, R. J., & Suri, J. S. (2013). Automated EEG analysis of epilepsy: A review. *Knowledge-Based Systems*, 45, 147-165.
- [11] Stam, C. J. (2005). Nonlinear dynamical analysis of EEG and MEG: Review of an emerging field. *Clinical Neurophysiology*, 116(10), 2266-2301.

- [12] Hazarika, N., Chen, J. Z., Tsoi, A. C., & Sergejew, A. (1997). Classification of EEG signals using the wavelet transform. *Signal Processing*, 59(1), 61-72.
- [13] Pigorini, A., Casali, A. G., Casarotto, S., Ferrarelli, F., Baselli, G., Mariotti, M., & Massimini, M. (2011). Time-frequency spectral analysis of TMS-evoked EEG oscillations by means of Hilbert–Huang transform. *Journal of Neuroscience Methods*, 198(2), 236-245.
- [14] Petrosian, A., Prokhorov, D., Homan, R., Dasheiff, R., & Wunsch, D. (2000). Recurrent neural network-based prediction of epileptic seizures in intra- and extracranial EEG. *Neurocomputing*, 30(1-4), 201-218.
- [15] Subasi, A., & Erçelebi, E. (2005). Classification of EEG signals using neural network and logistic regression. *Computer Methods and Programs in Biomedicine*, 78(2), 87-99.
- [16] Übeyli, E. D. (2009). Analysis of EEG signals by implementing eigenvector methods/recurrent neural networks. *Digital Signal Processing*, 19(1), 134-143.
- [17] Schirrmeister, R. T., Springenberg, J. T., Fiederer, L. D. J., Glasstetter, M., Eggensperger, K., Tangermann, M., ... & Ball, T. (2017). Deep learning with convolutional neural networks for EEG decoding and visualization. *Human Brain Mapping*, 38(11), 5391-5420.
- [18] Hosseini, M. P., Soltanian-Zadeh, H., Elisevich, K., & Pompili, D. (2016). Cloud-based deep learning of big EEG data for epileptic seizure prediction. *2016 IEEE Global Conference on Signal and Information Processing (GlobalSIP)*, 1151-1155.